

Welch Allyn RetinaVue™ 100 Imager



Directions for use

Software version 6.XX

WelchAllyn®

Advancing Frontline Care™

© 2018 Welch Allyn. All rights are reserved. To support the intended use of the product described in this publication, the purchaser of the product is permitted to copy this publication, for internal distribution only, from the media provided by Welch Allyn. No other use, reproduction, or distribution of this publication, or any part of it, is permitted without written permission from Welch Allyn. Welch Allyn assumes no responsibility for any injury to anyone, or for any illegal or improper use of the product, that may result from failure to use this product in accordance with the instructions, cautions, warnings, or statement of intended use published in this manual.

Software in this product is Copyright 2018 Welch Allyn or its vendors. All rights are reserved. The software is protected by United States of America copyright laws and international treaty provisions applicable worldwide. Under such laws, the licensee is entitled to use the copy of the software incorporated with this instrument as intended in the operation of the product in which it is embedded. The software may not be copied, decompiled, reverse-engineered, disassembled, or otherwise reduced to human-perceivable form. This is not a sale of the software or any copy of the software; all right, title, and ownership of the software remain with Welch Allyn or its vendors.

For information about any Welch Allyn product, contact Welch Allyn Technical Support: <http://www.welchallyn.com/support>.




411492, (printed copy)
DIR 80019780 Ver. J
Revision date: 2018-07

This manual applies to the  901103 Fundus Camera

Distributed by Welch Allyn, Inc.



Medimaging Integrated Solution Inc. (MiiS)
1F, No. 7, R&D Rd II, Hsinchu Science Park,
Hsinchu, TAIWAN 30076, R.O.C.

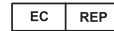

Regulatory Affairs Representative
MedNet GmbH
Borkstrasse 10
48163 Muenster
Germany

Made in Taiwan



WelchAllyn[®]
Advancing Frontline Care[™]

Contents

Symbols	1
Introduction	5
Intended use	5
About warnings and cautions	5
General warnings, cautions, and notes	5
Controls and connectors	8
Using the RetinaVue 100 Imager	11
Unpack the system	11
Docking station	17
Startup	18
First time selection of USB connectivity workflow	19
First time selection of Wi-Fi connectivity workflow	21
Select or change a clinic on the RetinaVue 100 Imager camera (wireless only)	28
Select or change a Multi-State clinic on the RetinaVue 100 Imager camera (wireless only)	29
Schedule a patient from the RetinaVue Network Customer Portal (wireless only)	31
Create a new patient on the RetinaVue 100 Imager camera (wireless and USB)	37
Select a patient from the patient list (wireless and USB)	40
Review patient tests (USB workflow)	42
About the Auto and Manual Exam modes	45
Adaptive software in the Auto exam mode for software 6.10.02 and higher	47
Perform an eye exam using the Auto exam mode	48
Perform a manual image capture while in Auto exam mode	54
Skip an image capture while in Auto exam mode	56
Perform a manual image capture while in Manual exam mode	57
Troubleshooting	59
Maintenance	69
Clean and disinfect	69
Specifications	71
Device radio	73
General compliance and standards	75
General radio compliance	75

Appendix	83
RetinaVue 100 Imager Home screen	83
Update the RetinaVue 100 Imager software to version 6.XX and above (USB and wireless)	97
Viewing or changing the RetinaVue 100 Imager settings	100
RetinaVue 100 Imager Accessories	106
Warranty	109

Symbols

Documentation symbols

For information on the origin of these symbols, visit <http://www.welchallyn.com/symbolsglossary> for the Welch Allyn symbols glossary.



WARNING The warning statements in this manual identify conditions or practices that could lead to illness, injury, or death. Warning symbols will appear with a grey background in a black and white document.



Caution The caution statements in this manual identify conditions or practices that could result in damage to the equipment or other property, or loss of data.



Consult directions for use (DFU). A copy of the DFU is available on this website. A printed copy of the DFU can be ordered from Welch Allyn for delivery within 7 calendar days.



Meets essential requirements of the European Medical Device Directive 93/42/EC.

Power symbols



Power button



Direct current (DC)



Rated power input, DC



Rated power output, DC



Fuse






Rechargeable battery

Li-ion








Battery Charging - Non-AC powered









Connectivity symbols

	USB		Wireless
	Expected connection to the the RetinaVue Network Server is not occurring.		

Wireless radio symbols

	Non-ionizing electromagnetic radiation		Wireless signal strength <ul style="list-style-type: none"> • Best (4 bars) • Good (3 bars) • Fair (2 bars) • Weak (1 bar) • No signal (no bars) • No connection (X)
	The identification number assigned by the Federal Communication Commission 2AFB3WA-WEC100		Industry Canada identification number. The equivalent governing body to the FCC in the United States 21375-WAWEC100 CAN ICES-3 (B)/NMB-3(B)
	Australian Communications and Media Authority (ACMA) Radio Compliance Mark (RCM)		

Shipping, storing, and environment symbols

	Temperature limit		Keep Dry
	Separate collection of Electrical and Electronic Equipment. Do not dispose as unsorted municipal waste.		Fragile
	Humidity limitation		Atmospheric pressure limitation
	This way up		Stacking limit by number



Recyclable

Miscellaneous symbols



Product Identifier



Reorder Number



Global Trade Item Number



For indoor use only



Type BF applied part




Prescription only or "For Use by or on the order of a licensed medical professional"



Clock; time switch; timer



Serial Number

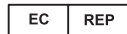

Authorized Representative in the European Community



Degree of protection provided by enclosure (IP code)

Battery charge level indicators



Device is running on battery and charge is between 75 and 100 percent.



Device is running on battery and charge is below 75 percent.



Device is running on battery and charge is below 50 percent.



Device is running on battery and charge is below 25 percent.



Device is running on battery and charge is below 10 percent. Device blinks with alternating icons.



Device is connected to a power supply and is docked.

Introduction

Intended use

The Welch Allyn RetinaVue 100 Imager is a digital hand-held eye-fundus camera used to record digital photographs of the fundus (including retina, macula and optic disc) of the human eye and surrounding area.

Prescription only or "For Use by or on the order of a licensed medical professional".

About warnings and cautions

Warning and caution statements can appear on the RetinaVue 100 Imager device, the packaging, the shipping container, or in this *Directions for use*.

The RetinaVue 100 Imager is safe for patients and clinicians when used in accordance with the instructions and the warning and caution statements presented in this *Directions for use*.

Before using the device, you must familiarize yourself with all warnings and cautions, with the steps to power up the device, and with the sections of this *Directions for use* that pertain to your use of the device. Specific warnings and cautions are also found throughout this manual.

- Failure to understand and observe any warning statement in this manual could lead to patient injury or illness.
- Failure to understand and observe any caution statement in this manual could lead to damage to the equipment or other property, or loss of patient data.

General warnings, cautions, and notes



WARNING Exposure to the RetinaVue 100 Imager flash, in rare cases, may contribute to migraines or temporary visual disturbances, particularly in patients with a history of light sensitivity or migraines. Multiple flashes per eye in a single exam may increase the likelihood of occurrence.



WARNING Clean and disinfect the eye cup after each patient to avoid the risk of cross-contamination.



WARNING Operating the device in the presence of a flammable anesthetic mixture with air, oxygen, or nitrous oxide; in oxygen-enriched environments; or in any other potentially explosive environment may cause a fire and explosion hazard.



WARNING Do not use batteries from other sources. This may result in excessive temperatures, fire, or explosion.



WARNING To prevent fire or electrical shock, do not expose the device to rain or moisture.



CAUTION When manually entering patient information, do not use the same patient identifier for multiple patients. This may result in inaccurate data association.



CAUTION Avoid subjecting the device to vibration or shock.



CAUTION Avoid using the device in a dusty environment.



CAUTION Keep the power cord away from any heat source.



CAUTION Do not sterilize the device.



CAUTION All servicing of this equipment is to be performed by Welch Allyn, Inc. or by a center approved by Welch Allyn, Inc.



CAUTION This device needs to be installed and put into service according to the EMC information provided herein.



CAUTION Portable and mobile RF communications equipment can affect the RetinaVue 100 Imager performance.



CAUTION Input overload can occur in proximity to defibrillator or electrocautery devices.



CAUTION This device is not protected against the ingress of water and should not be used in the presence of liquids which may enter the device.



CAUTION Results deleted on the device cannot be recovered.



CAUTION Do not modify this equipment without authorization of the manufacturer.



CAUTION Do not use the device in direct sunlight. Strong ambient light may affect results.



CAUTION Within the patient environment, only connect the equipment (the RetinaVue 100 Imager) to computers or laptops that have passed the safety standard for information technology equipment IEC 60601-1, EN 60601-1, ANSI/AAMI ES 60601-1 to ensure the safety of the USB electrical connection. Outside the patient environment, only connect the equipment (the RetinaVue 100 Imager) to computers or laptops that have passed the safety standard for information technology equipment IEC 60950-1, EN 60950-1, UL 60950-1 to ensure the safety of the USB electrical connection.



CAUTION Do not turn off the device while saving data from a patient test. Turning off the device while saving patient data may delete the data.



CAUTION Do not use the device adjacent to or stacked with other equipment. If the device is used adjacent or stacked with other equipment, the device should be observed to verify normal operation in the configuration in which it will be used.



CAUTION Do not submerge the device in liquid or allow liquid to enter the interior of the device.



CAUTION Do not use automatic cleansing machines or sterilization.



CAUTION Consult the cleaning agent and germicidal cleaner agent manufacturer instructions for their proper use and germicidal efficacy prior to their use.




CAUTION Only use the cleaning or germicidal cleaner agent types listed or damage may occur.

Note

The RetinaVue 100 Imager can be used on patients with eye implants and intraocular lenses. While the RetinaVue 100 Imager can be used on patients with cataracts and other eye opacities, the use of the RetinaVue 100 Imager may result in a lower quality image due to the increased reflection of the flash off the patient's intraocular lens.

Controls and connectors

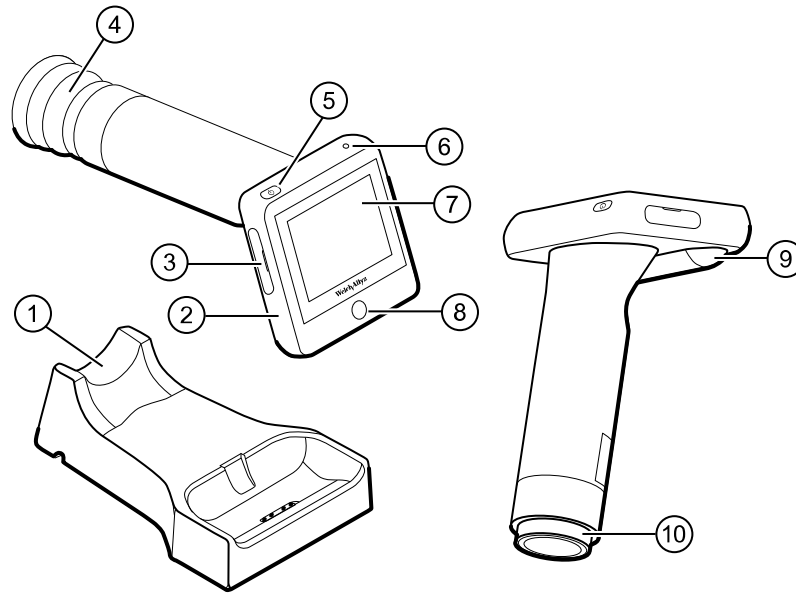

No.	Feature	Description
1	Docking station	<p>Charges the RetinaVue 100 Imager Pro and enables data transfer to a computer. Connect to an electrical outlet or to the “Always On” USB port of a computer.</p> <p>Note USB ports supply power only when the computer is on. However, many newer computers contain a yellow “Always On” USB port that can supply power even when the computer is sleeping, hibernating, or turned off. See you computer manufacturer’s user documentation for instructions on changing the power management options and for identifying the “Always On” USB port.</p>
2	RetinaVue 100 Imager device	Digital hand-held eye-fundus camera used to record digital photographs of the fundus of the human eye and surrounding area.
3	SD Card port/service port	<p>SD Card slot. Insert SD Card for data storage.</p> <p>Note The service port under the SD Card cover is for factory use only and is not for customer use.</p>
4	Eye cup	Blocks ambient light during testing
5	Power button	<ul style="list-style-type: none"> • Push button to turn the device on and off. • Push button to bring the device out of Sleep mode. • Push button to launch a dialogue box with three options: <ul style="list-style-type: none"> ○ Power off ○ Sleep


No.	Feature	Description
		○ Cancel
6	Power indicator	<ul style="list-style-type: none">• When the LED light is green, it indicates that the device is docked and is fully powered.• When the LED light is amber, it indicates that the device is docked and is in the process of charging.• When the device is not docked the LED is off.
7	LCD touchscreen	3.5" color HD LCD touchscreen. Displays information for charging, testing, and adjustments to the device.
8	Home button	Press the Home button to navigate to the Home screen (workflow selection) from any screen on the device.
9	Battery cover	Covers the battery compartment
10	Bezel	Area to which the eye cup attaches

Using the RetinaVue 100 Imager

View images saved to the RetinaVue 100 Imager on the LCD touchscreen. The images saved on the SD card are encrypted. To view the images on the PC, install the RetinaVue Network software. For information on the RetinaVue Network Software visit www.RetinaVue.com.

Unpack the system

The RetinaVue 100 Imager device is packaged with the following items. Confirm that all items are included.

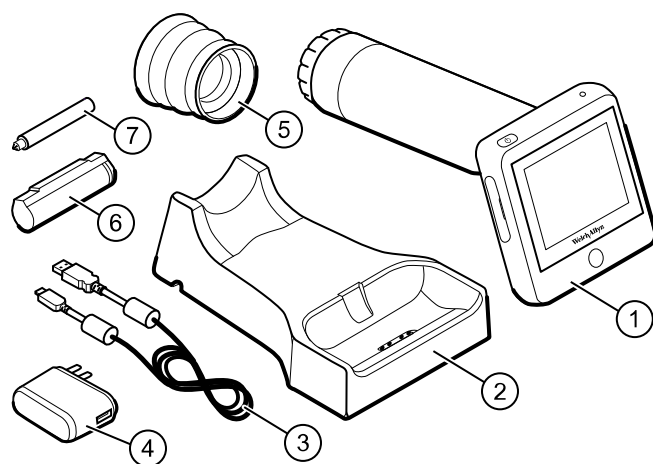

Number	Item	Description
1	RetinaVue 100 Imager device	Digital hand-held eye-fundus camera used to record digital photographs of the fundus of the human eye and surrounding area.
2	Docking station	Charges the RetinaVue 100 Imager device and enables USB data transfer to a computer or laptop. Connect to region-specific AC wall outlet.
3	USB cable	Connects the docking station to a computer or laptop to transfer results. Connects the docking station to the wall adapter plug for AC charging.
4	Wall adapter plug	Connects to USB cable. Match to region-specific AC wall outlet.

Number	Item	Description
5	Eye cup	Blocks ambient light during testing
6	Battery	Rechargeable lithium ion battery
7	Screwdriver	Phillips screwdriver to remove battery cover screws
Not shown	RetinaVue 100 Imager <i>Directions for use</i>	This document. The <i>Directions for use</i> (DFU) is also available on the Welch Allyn web site.
Not shown	SD Card	SD Card installed in the SD Card slot for data storage
Not shown	Lens spray	RetinaVue 100 Imager lens cleaner
Not shown	Lens cloth	RetinaVue 100 Imager lens cleaning cloth

If any item is missing, contact Welch Allyn Technical Support: www.welchallyn.com/service. We recommend that you save the shipping box and packing materials in case you need to store or ship the system.

Battery installation and replacement

Use the device on battery power once the battery is installed and charged. The battery is charged through contact pins on the device and the docking station.

Install the battery/replace the battery

Set the RetinaVue 100 Imager on a flat, clean surface with the LED display facing down.




WARNING Personal injury risk. Improper handling of the battery can lead to heat generation, smoke, explosion, or fire. Do not short-circuit, crush, incinerate, or disassemble the battery. Never dispose of batteries in refuse containers. Always recycle batteries according to national or local regulations.




WARNING Use only Welch Allyn approved accessories, and use them according to the manufacturer's directions for use. Using unapproved accessories with the camera can affect patient and operator safety and can compromise product performance and accuracy, and void the product warranty.


1. Loosen the screw at the left side of the battery cover.
2. Loosen the screw at the right side of the battery cover.
3. Push in the battery cover tab and then lift up the battery cover. Slide open the battery cover by tilting the cover upward.

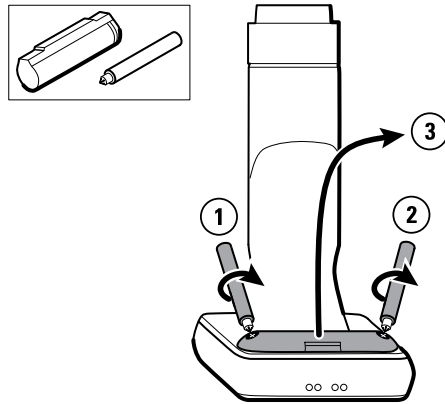

4. For first-time setup of the RetinaVue 100 Imager, insert the battery into the battery compartment. For battery replacement, remove the existing battery and then insert the new replacement battery into the battery compartment.
 - a. This procedure applies to first-time setup of the RetinaVue 100 Imager. Insert the battery into the battery compartment. Match the polarity of the battery to the + and - markings on the RetinaVue 100 Imager.
 - b. This procedure only applies to battery replacement. Remove the existing battery from the battery compartment and insert the new replacement battery into the battery compartment. Match the polarity of the battery to the + and - markings on the RetinaVue 100 Imager.
5. Rotate the raised flat edge of the battery until it sets in the back notch of the battery compartment.

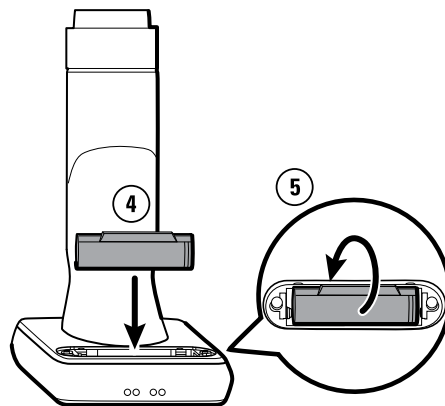

6. Insert the two battery cover tabs into the slot, then push in the tab to close the battery compartment.
7. Tighten the screw at the left side of the battery cover.
8. Tighten the screw at the right side of the battery cover.

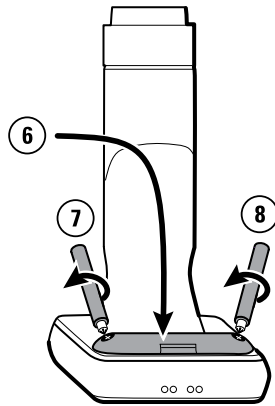


Note Do not over-tighten the screws.

Eye cup


The eye cup blocks ambient light during testing and covers the patient's eye socket to allow them to best focus on the green fixation lights. The eye cup works on the patient's left and right eye.

Set up the eye cup for a patient exam

1. Remove the lens cap from the RetinaVue 100 Imager.

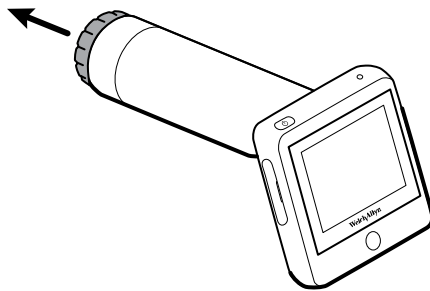

2. Slide the eye cup over the bezel on the RetinaVue 100 Imager.

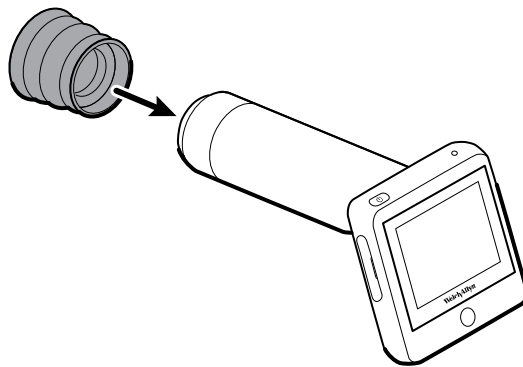

3. Use the eye cup to cover the patient's examined eye and allow them to focus on the green fixation lights.



Hold the RetinaVue 100 Imager for a patient exam

The RetinaVue 100 Imager must be held level and in line with your patient's eye.

There are two recommended techniques to properly hold the RetinaVue 100 Imager for a quick and successful retinal exam:

Telescope Grip



Smart phone Grip

Make micro adjustments with your guide hand to fine tune the focus of the RetinaVue 100 Imager.

Guide (Patient End) Hand


Hold the patient end of the RetinaVue 100 Imager with just the thumb and forefinger. Use the small and ring finger to steady the RetinaVue 100 Imager by anchoring them gently on the patient's forehead.

Docking station


Place the device into the docking station to charge the device or to connect and transfer data to a computer or laptop PC running the RetinaVue Network™ software. (Requires additional software).

Set up the docking station for charging


1. Turn the docking station over.

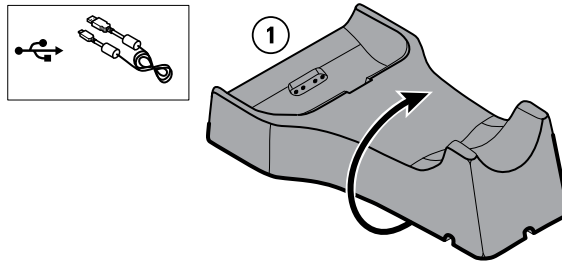


2. Connect the USB cable to the docking station.
3. Either plug the USB cable into the wall adapter plug and then plug the wall adapter into an AC outlet, or plug the USB cable into the "Always On" USB port on a computer.
 - a. Plug the USB cable into the wall adapter plug and then plug the wall adapter into an AC outlet.

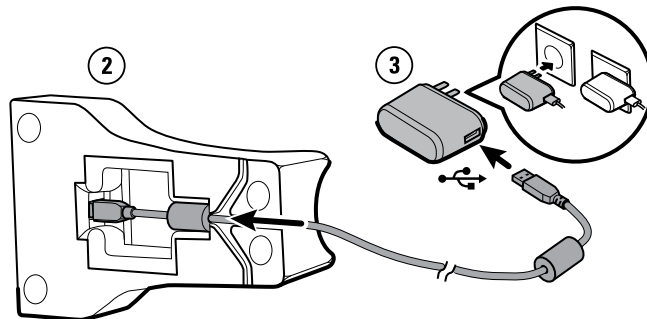

- b. (Optional) Plug the USB cable into the "Always On" USB port on a computer.

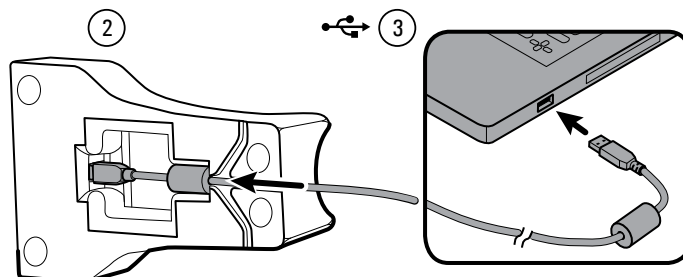


4. Place the RetinaVue 100 Imager on the docking station charger with the contact pins of the RetinaVue 100 Imager touching the contact pins of the docking station charger.

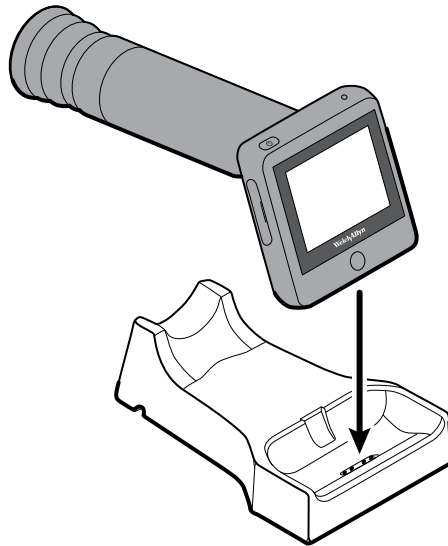

5. Before the first use, allow the RetinaVue 100 Imager to charge for 5 hours on the docking station charger.

Note The power indicator light turns amber to indicate that the RetinaVue 100 Imager is in the processing of charging. If the power indicator light does not turn on, ensure that the contact pins of the RetinaVue 100 Imager are touching the contact pins of the docking station charger and reseal the RetinaVue 100 Imager into the docking station. Charging is impacted if the pins are not properly touching.

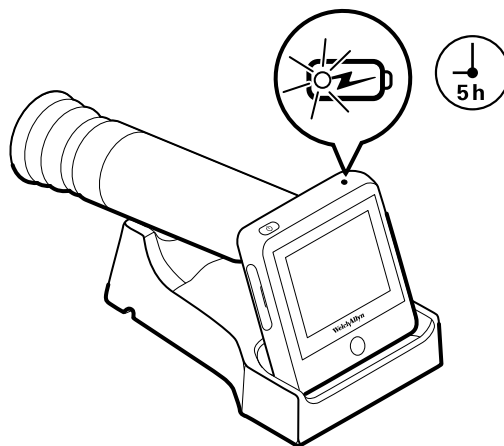


Note The removal of the plug from the socket provides the means of isolation from the mains power supply. The user is reminded to keep this means of disconnection easily accessible in case of emergency.

Startup

Power

The Power button, located on the upper-left corner of the RetinaVue 100 Imager, performs multiple functions:

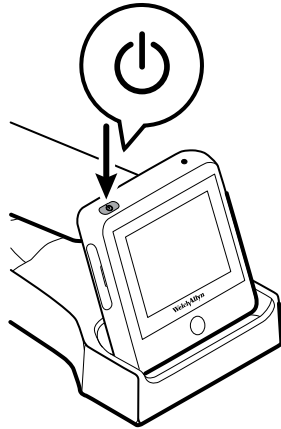

- powers up the RetinaVue 100 Imager
- wakes the RetinaVue 100 Imager from Sleep mode
- opens a dialogue screen with controls to power off, enter Sleep mode, or to cancel and return to the prior screen, depending on the option you choose


Power up the RetinaVue 100 Imager

The RetinaVue 100 Imager runs a brief diagnostic self-test each time it powers up.



CAUTION Always use the RetinaVue 100 Imager with an adequately charged and properly functioning battery.

Press  to power up the RetinaVue 100 Imager.

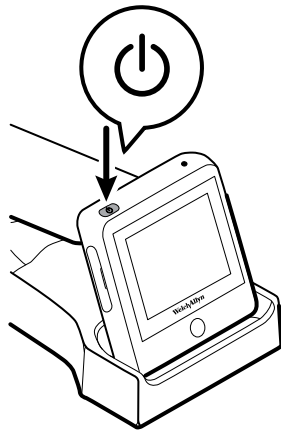


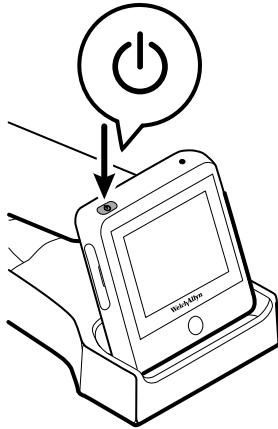
As the device is powering up, the RetinaVue 100 Imager displays the Welch Allyn logo.

- If your facility has chosen to use a Device access code or an Advanced settings code, then the first image you see is the login screen.
- If your facility has not chosen to use a Device access code or Advanced settings code, then the first image you see is the Home screen.

First time selection of USB connectivity workflow

The RetinaVue 100 Imager runs a brief diagnostic self-test each time it powers up.

1. Press  to power up the RetinaVue 100 Imager.



After the RetinaVue 100 Imager displays the Welch Allyn logo, the first-time camera setup wizard screen appears.

Note The USB or Wi-Fi camera setup only needs to be completed once unless a change to the connectivity workflow is desired.

2. Touch **RetinaVue Network - USB** to select the USB workflow and touch **Next**.




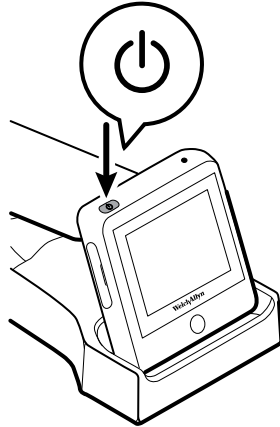
Note For the USB workflow, the RetinaVue 100 Imager must be connected to a computer running the RetinaVue Network software (requires additional software) to transfer exams.

Note To change the connectivity workflow settings after the initial set up, touch **Settings > Advanced settings > Restore factory defaults**. Touch **OK** at the dialog: *Reset all the settings to defaults?* and then touch **OK** at the dialog *All settings were restored to factory defaults.*

First time selection of Wi-Fi connectivity workflow

The RetinaVue 100 Imager runs a brief diagnostic self-test each time it powers up.

1. Press  to power up the RetinaVue 100 Imager.







After the RetinaVue 100 Imager displays the Welch Allyn logo, the first-time camera setup wizard screen appears.

Note The USB or Wi-Fi camera setup only needs to be completed once unless a change to the connectivity workflow is desired.

2. Touch **RetinaVue Network - Wi-Fi** to select the wireless workflow and touch **Next**.



3. Touch  (up) or  (down) to enter, or adjust, the Day, Month, and Year. Touch **Next**.
4. Touch  (up) or  (down) to enter, or adjust, the Hour, Minute (Min), Continent, and Location. Touch **Next**.

Note Use the radio buttons to select the 24-hour or 12-hour format. The available wireless networks appear on the *Select a Wi-Fi connection* screen.

5. From the list of available networks, touch on the Wi-Fi connection that you want to connect to.



Note *Wi-Fi network names shown here are examples.

Note If you don't see the desired Wi-Fi network that you want to connect to, you may need to scroll down.


Note During Wi-Fi communication, the RetinaVue 100 Imager camera only supports the use of the following secure authentication types: WPA2, WPA2-PSK, WPA2-Enterprise (EAP-PEAP), and IEEE 802.11i. Open networks from public access points are not supported. See the *Device radio* section for further information. If you don't see the network name in your list it might be a hidden network. Select "Other" and enter the network name. See the *Troubleshooting* section for further information.

6. Using a secured network, enter your network Wi-Fi Username and Password or a Passphrase and touch **Connect**.
 - a. (Optional). Touch **Remember this connection** and touch **Connect** to use the selected network for subsequent Wi-Fi connections.
7. To continue setup, launch an internet browser from a computer and enter the following address: https://www.retinavue.net/RN_CustomerPortal/.
8. On your RetinaVue 100 Imager camera, touch **Next**.
9. On the RetinaVue Network Customer Portal, enter your RetinaVue Network Username and Password then click **Log In**.



Note If prompted, click **I agree** to accept the End User License Agreement on the RetinaVue Network Customer Portal.

- On the RetinaVue Network Customer Portal, select *Manage Devices (wireless only)* and click **Next**.

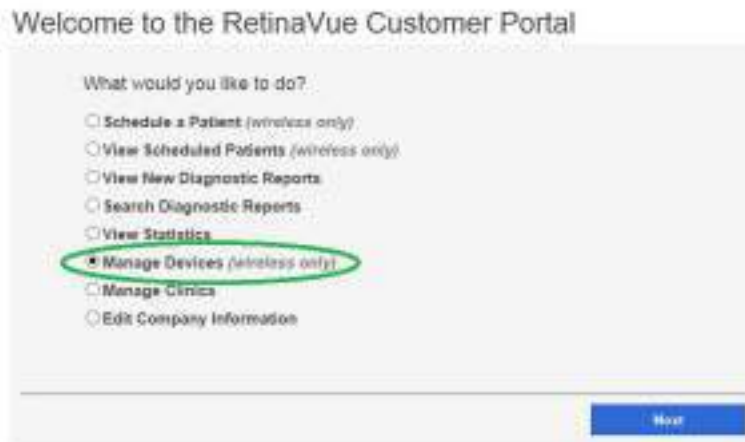

- Click **Add Device**.



- On your RetinaVue 100 Imager camera, touch **Next** to generate a Device Registration Code.
- On the RetinaVue Network Customer Portal, enter the Device Registration Code from the RetinaVue 100 Imager and then click **Enter**.



Note *The Device Registration Code shown here is an example.

Note The *Next* button is greyed out until the Device Access Code is received by the RetinaVue 100 Imager camera from the RetinaVue Network Customer Portal.

Note The Device Registration Code needs to be entered into the RetinaVue Network Customer Portal as it appears on the RetinaVue 100 Imager *device installation* screen including the dash (-).

Note The code expires after 30 minutes.

14. On the RetinaVue Network Customer Portal , enter a Device Name (optional) and select a clinic, or clinics, where the RetinaVue 100 Imager camera will to be used. Select each clinic by clicking the check box next to the clinic name and then click **Add Device**.

Warning: This registration request will expire if the device is not added by 10/3/2017 10:05:31 AM -04:00

Add New Device

Device Name
Mountain View And Valley View Clinic

Device Serial Number
28102004

Device Access Code
####

Select clinics to associate with this camera


	Clinic Name	Address
<input checked="" type="checkbox"/>	Mountain View Clinic	1234 Red River Road, Canyon Top, UT
<input checked="" type="checkbox"/>	Valley View Clinic	12345 Riverbent Ave, Valley View, UT

Back Add Device


- (Optional) To select multiple clinics where the RetinaVue 100 Imager camera will be used at, select each clinic by clicking the check box next to the clinic name.
- (Optional) To enable a Device Access Code for the RetinaVue 100 Imager camera, select the **Device Access Code** checkbox and enter your code into the *Device Access Code* field. The camera and the RetinaVue Network Customer Portal support from 4 to 6-digits in the Device Access Code.

Note Setting the Device Access Code before first use provides additional security features in the RetinaVue 100 Imager to prevent unintended access to patient information by unauthorized individuals.


- On your RetinaVue 100 Imager camera, touch **Next** at the *Device Registration Code* screen.
- On your RetinaVue 100 Imager camera, touch **Next** on the screen: *Camera registered successfully*.

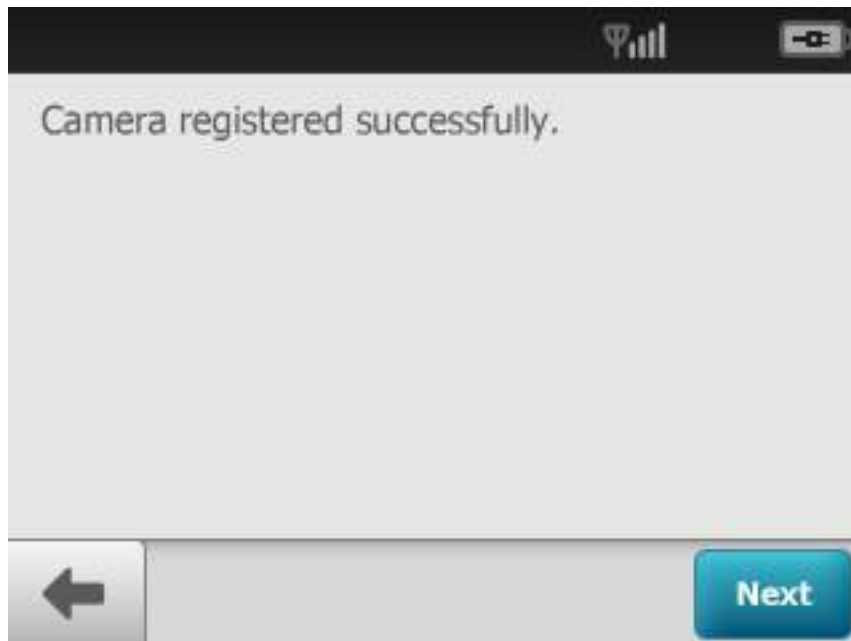


17. (Optional) If a Device Access Code was selected at the RetinaVue Network Customer Portal, touch **OK** at the dialogue: *A new Device Access Code has been downloaded. The camera will require this code on next encounter.*

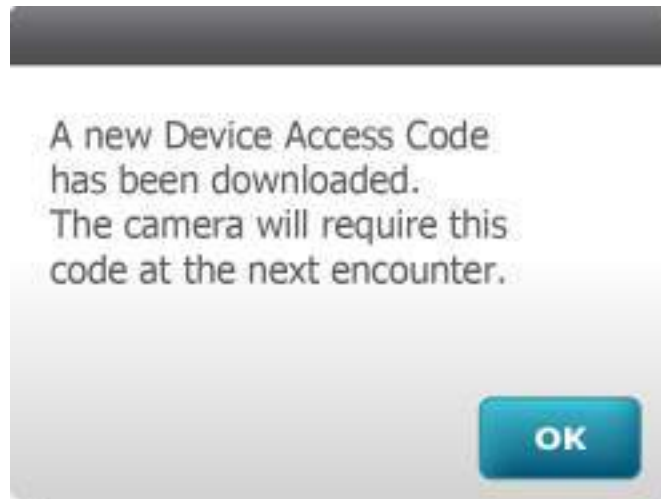

18. (Optional) If multiple clinics were selected at the RetinaVue Network Customer Portal, use the drop-down menu to select the primary clinic where exams will be performed and then touch **Next**.

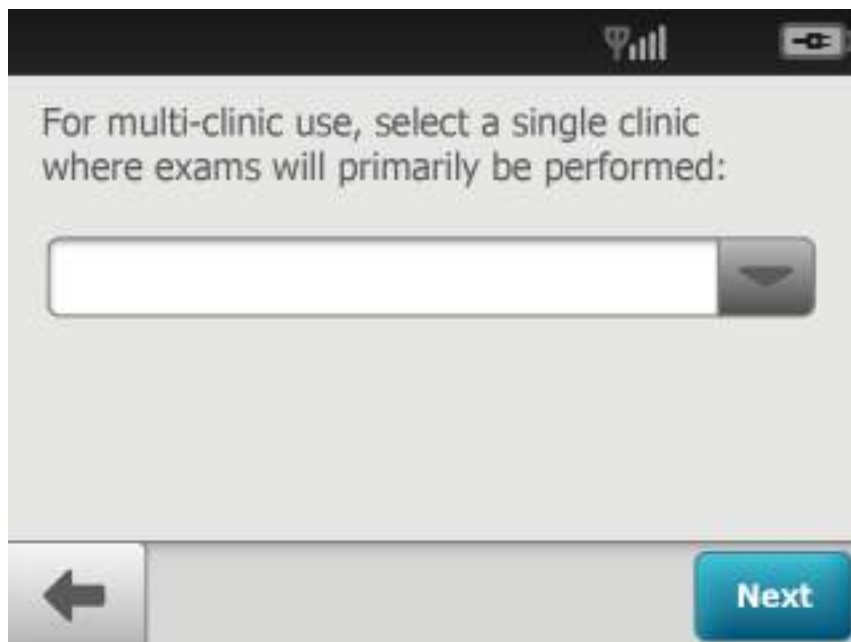

19. Touch **OK** at the *Setup complete* screen.

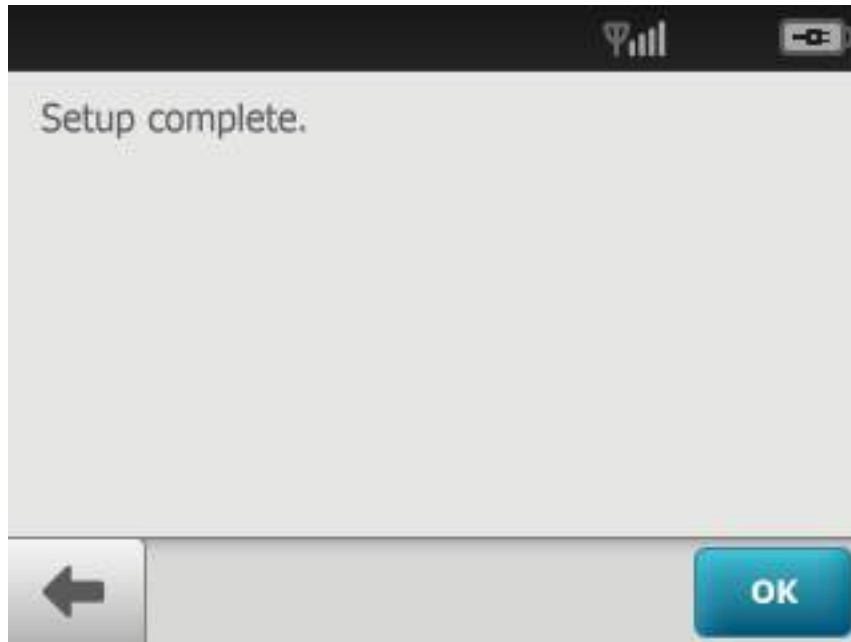



20. (Optional) If a Device Access Code was set at the RetinaVue Network Customer Portal, enter the Device Access Code on the RetinaVue 100 Imager and touch **OK**.



Select or change a clinic on the RetinaVue 100 Imager camera (wireless only)


1. Touch **Settings** on the RetinaVue 100 Imager *Home* screen.
2. Touch **Clinic / States**.
3. Use the drop-down menu to select a clinic.



4. Touch  (Back) to return to the *Settings* screen.
5. Touch  (Home) to return to the RetinaVue 100 Imager *Home* screen.



Select or change a Multi-State clinic on the RetinaVue 100 Imager camera (wireless only)

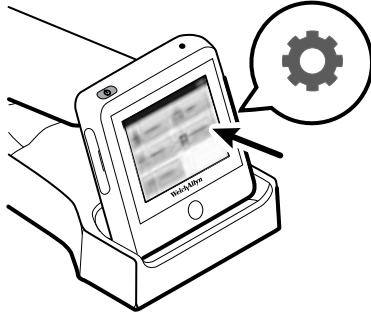1. Touch **Settings** on the RetinaVue 100 Imager Home screen.



2. Touch **Clinic / States**.
3. If multiple clinics are set up, use the drop-down menu to select the clinic you would like to set up for Multi-State use.
4. Touch **Multi-State**. The switch changes from Off to On.




5. Touch **States** to select between 2 and 6 States for the clinic. Touch **More** to see additional States and to proceed to OK when the desired number of States are selected.
6. Touch **OK** to return to the *Clinic/State Settings* screen.
7. Touch  (Back) to return to the *Settings* screen.
8. Touch  (Home) to return to the RetinaVue 100 Imager *Home* screen.

Note For each exam performed in multi-state use, select the state where the exam was performed. Touch **OK** to confirm.



Schedule a patient from the RetinaVue Network Customer Portal (wireless only)

1. Select **Schedule a Patient (wireless only)** from the RetinaVue Network Customer Portal *Welcome screen* and then click **Next**.

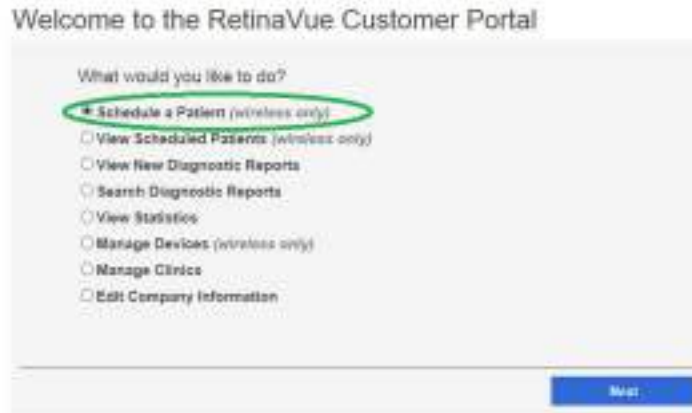


2. Use the drop-down menu to select a clinic.



The *Schedule a Patient* screen appears for your clinic.

3. Complete the required fields and click **Schedule Patient**.

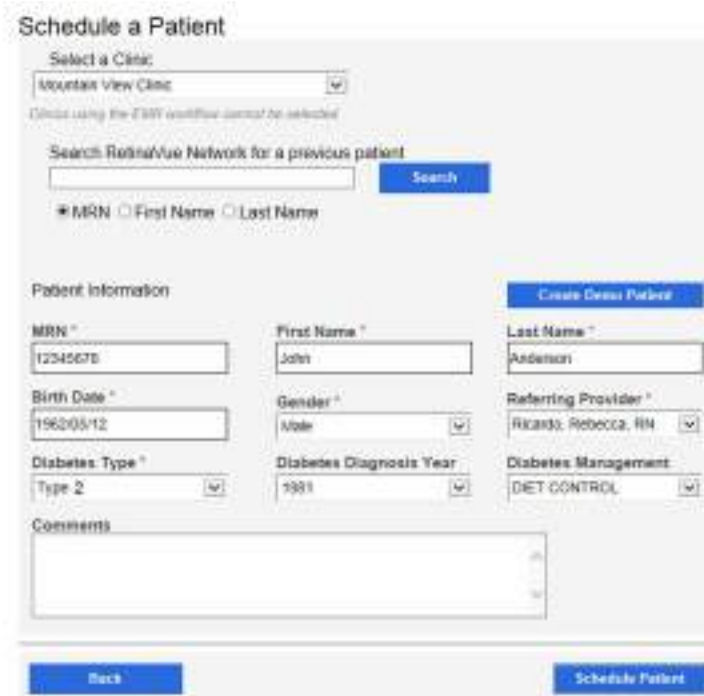




CAUTION When manually entering patient information, do not use the same patient identifier for multiple patients. This may result in inaccurate data association.


Note Required fields contain a red asterisk. Required fields include: MRN, First and Last Name, Date of Birth, Gender, Referring Provider, Diabetes Type. Optional fields include: Diabetes Diagnosed and Diabetes Management.

Note To search for a patient who has had images captured in the past, select MRN, First Name, or Last Name and enter either the patient's MRN, or their First Name, or their Last Name in the Search field and click **Search**. Any previous patients of your clinic (within the RetinaVue Network) meeting any of the search parameters appear on the *Patients Search* screen.


4. Ensure that the RetinaVue 100 Imager camera is powered on and connected to a secure Wi-Fi network. Touch **Patients**. The *Patients screen* appears.



- (Optional) Touch  to search for a patient. Select a patient. Touch on a patient's name to select that patient. The *Patients* screen appears.

- (Optional) Touch  to refresh the patient list from the RetinaVue Network Customer Portal Scheduled Patients list.

5. Confirm the patient demographics and touch **More** to view the second *Patient summary* screen with additional Referring Physician and Diabetes type required fields.

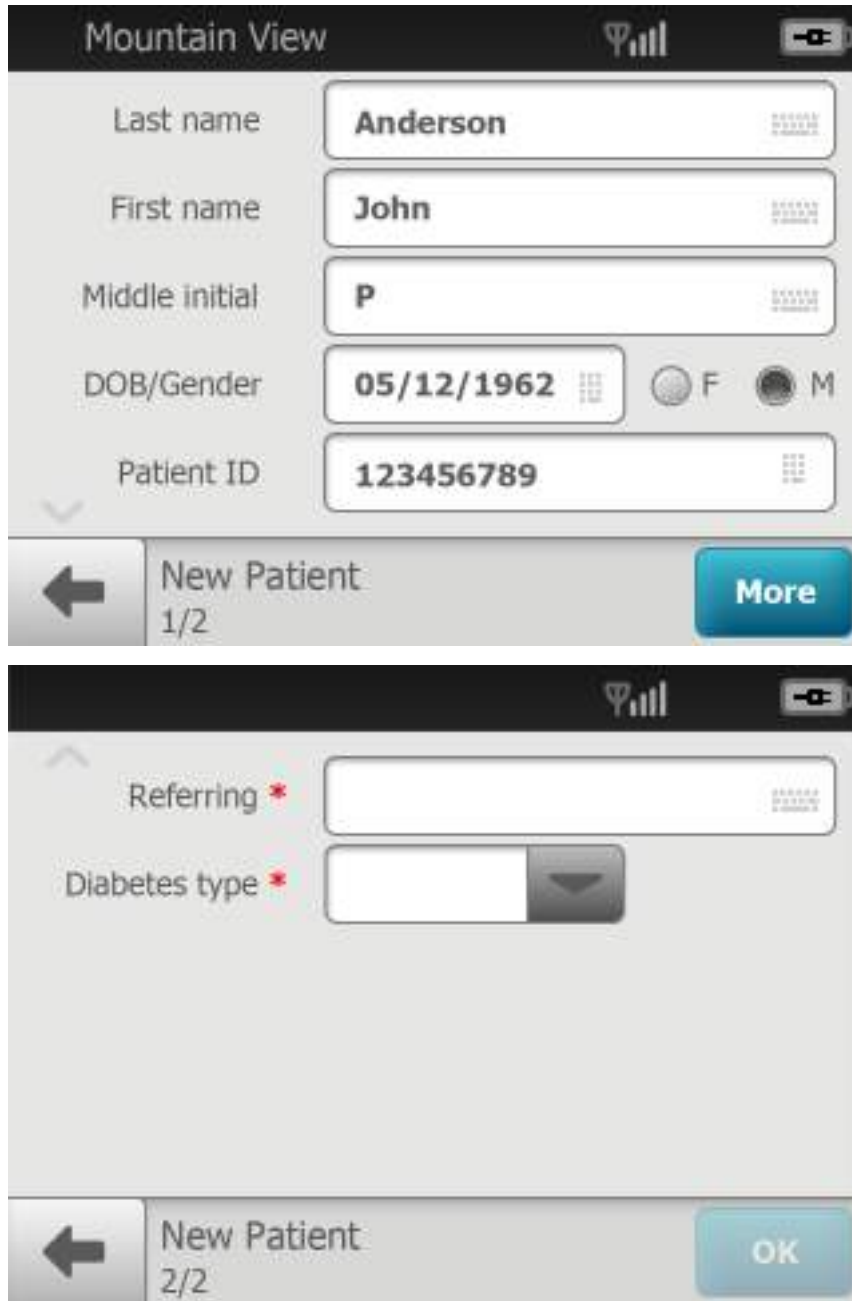


6. Touch **OK** to start a new exam. The *Pre-acquisition* screen appears.



7. Refer to the section *Perform an eye exam using the Auto exam mode* for instructions on performing exams and acquiring images.

View Scheduled Patients on the RetinaVue Network Customer Portal (wireless only)

1. Select **View Scheduled Patients (wireless only)** from the RetinaVue Network Customer Portal *Welcome screen* and then click **Next**.

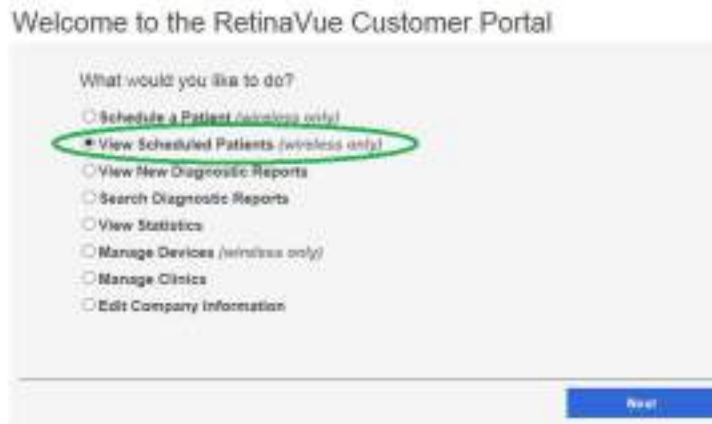

The *View Scheduled Patients* screen appears.

2. Use the drop-down menu to select a clinic and then view scheduled patients.

View Scheduled Patients

Mountbain View Clinic

[Schedule Patient](#) [View Reports](#)

Clicking on the ID number cannot be selected

Patient Name	MRN	Birth Date	Status	
Demo Patient	DCMD065T88773195	12/05/1973	Scheduled	
Demo Patient	DCMD065T88773196	12/05/1973	Submitted	
Jim Smith	12345678	8/20/2000	Scheduled	
Carole Winespinner	12345679	3/1/1984	Scheduled	
Wendy Cooper	12345680	8/17/1970	Scheduled	
John Anderson	12345681	5/12/1962	Scheduled	

[Back](#)

Note Scheduled patients appear as *Scheduled* in the Status column. If the patient's exam has been completed the patient record appears as *Submitted* until the patient's report is returned.

Schedule a Demo patient from the RetinaVue Network Customer Portal without submitting an exam for over-read (wireless only)

1. Select **Schedule a Patient (wireless only)** from the RetinaVue Network Customer Portal *Welcome* screen and then click **Next**.

Welcome to the RetinaVue Customer Portal

What would you like to do?

- Schedule a Patient (wireless only)**
- View Scheduled Patients (wireless only)
- View New Diagnostic Reports
- Search Diagnostic Reports
- View Statistics
- Manage Devices (wireless only)
- Manage Clinics
- Edit Company Information

[Next](#)

2. Use the drop-down menu to select a clinic.

Schedule a Patient

Select a Clinic:

[Back](#)

The *Schedule a Patient* screen appears for your clinic.

3. Click **Create Demo Patient**, complete all required fields, and then click **Schedule Patient**.

Schedule a Patient

Select a Clinic
Mountain View Clinic

Clinics using the ICDM workflow cannot be selected

Search RetinaVue Network for a previous patient

MRN First Name Last Name

Patient Information

MRN * First Name * Last Name *

Birth Date * Gender * Referring Provider *

Diabetes Type * Diabetes Diagnosis Year Diabetes Management

Comments

- The *View Scheduled Patients* screen appears.

View Scheduled Patients

Mountain View Clinic

* Clinics using the ICDM workflow cannot be selected

Patient Name	MRN	Birth Date	Status	
Demo Patient	DEM0394766775195	12/25/1225	Scheduled	X
Demo Patient	DEM0992301811835	12/25/1225	Scheduled	X
Jim Sentry	12345678	02/02/2000	Scheduled	X
Carroll Winesapper	12345678	3/1/1964	Scheduled	X
Woody Cuppans	12345678	9/17/1970	Scheduled	X
John Anderson	122915078	5/12/1952	Scheduled	X

- Ensure that the RetinaVue Imager camera is powered on and connected to a secure Wi-Fi network. Touch **Patients**. The *Patients* screen appears.




6. Touch on **Patient, Demo** to select the Demo Patient. The *Patient Summary* screen appears.
7. Confirm the demo patient demographics and touch **More** to view the second *Patient summary* screen with additional Referring Physician and Diabetes type required fields.
8. Touch **OK** to start a new exam to practice with a demo patient without submitting the exam for over-read on the RetinaVue Network Customer Portal.
9. Refer to the section *Perform an eye exam using the Auto exam mode* for instructions on performing exams and acquiring images.

Create a new patient on the RetinaVue 100 Imager camera (wireless and USB)


1. From the *Home screen*, touch **Patients**. The *Patients screen* appears. (**Wireless workflow images shown**)



(Optional for wireless) Touch  to search for a patient.

(Optional for wireless) Touch  to refresh the patient list from the RetinaVue Network Customer Portal Scheduled Patients list.

2. Touch  (New patient). The *New Patient* screen appears. Required fields are noted with a red asterisk.

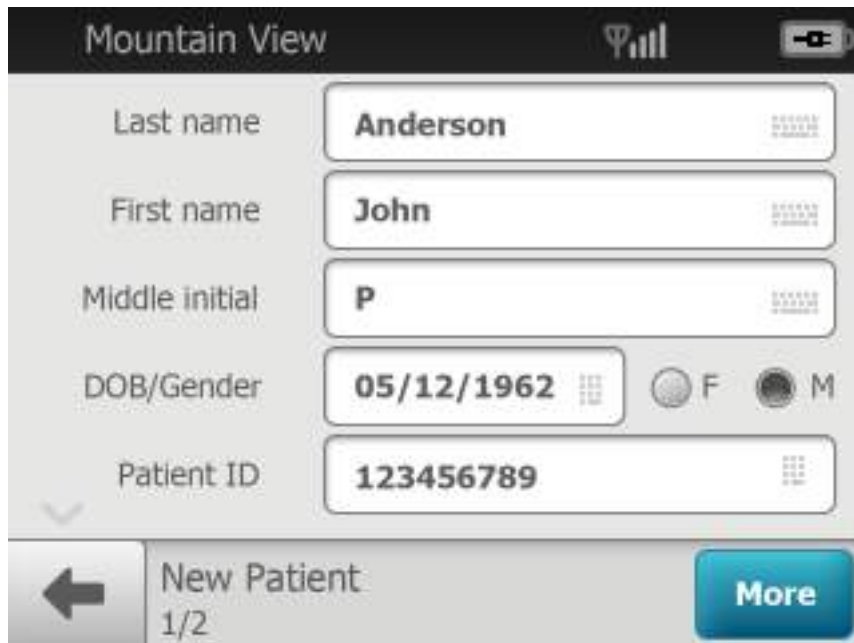



CAUTION When manually entering patient information, do not use the same patient identifier for multiple patients. This may result in inaccurate data association.

3. Touch within the Patient's Last name, First name, or Middle initial entry box to launch the on-screen keyboard.
4. Enter the patient's last name and first name. Touch **OK**. (Middle initial is optional.)
5. Enter the patient's Date of Birth. (Touch within the DOB/Gender field to launch the on-screen keyboard and to enter the patient's Date of Birth.) Touch **OK**.
6. Enter the patient's gender using the radio button to select male or female.
7. Enter the patient's ID. Touch **OK**.
8. (Wireless only) Touch **More** to view the second *New Patient* screen with additional Referring Physician and Diabetes type required fields.

The screenshot shows a mobile application interface for entering patient information. At the top, there are status icons for signal strength and battery. Below that, there is a 'Referring' field with the text 'Murphy MD, Ryan' and a small 'More' icon to its right. Below the 'Referring' field is a 'Diabetes type' field with a dropdown menu showing 'II'. At the bottom of the screen, there is a navigation bar with a back arrow on the left, the text 'Patient summary' in the center, and a blue 'OK' button on the right.

9. Once all the required fields are complete, touch **OK** to open a new exam. The *Pre-acquisition* screen appears.



10. Refer to the section *Perform an eye exam using the Auto exam mode* for instructions on performing exams and acquiring images.

Select a patient from the patient list (wireless and USB)

1. From from the *Home screen*, touch **Patients**. The *Patients screen* appears. **(Wireless workflow images shown)**



2. Touch the Patient's name to view a summary of the patient information on the *Patient Summary screen*.

The image displays two screenshots of the 'Mountain View' patient entry interface. The top screenshot shows the 'Patient summary' screen with the following fields: Last name (Anderson), First name (John), Middle initial (P), DOB/Gender (05/12/1962, M), and Patient ID (123456789). A 'More' button is visible. The bottom screenshot shows the same screen with additional fields for Referring (Murphy MD, Ryan) and Diabetes type (II), and an 'OK' button.

3. Confirm the patient demographics and touch **More** (wireless only) to view the second *Patient summary* screen with additional Referring Physician and Diabetes type required fields.
4. Touch **OK** to start a new exam. The *Pre-acquisition* screen appears.




5. Refer to the section *Perform an eye exam using the Auto exam mode* for instructions on performing exams and acquiring images.



Review patient tests (USB workflow)

1. From the *Home screen*, touch **Review**. The patient's *Review screen* appears.

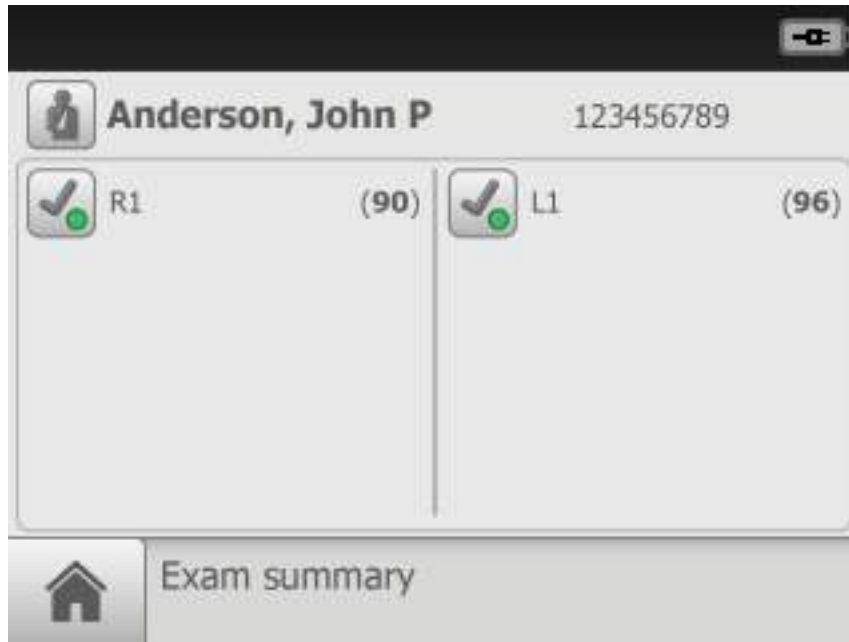

2. Touch on a patient's information to review the patient demographics and to view images associated with a patient.
3. Any images which have not been sent for overread appear on the *Exam summary screen*.



4. After the review is complete, touch  (Back) to return to the *Review screen* and then touch  (Home) to return to the *Home screen*.

About image quality scores

**Good quality
green**

If the image is of good quality, it is very likely to be read by a specialist.
(Scores of 40 +)



**Acceptable quality
yellow**

If the image is of good enough quality, it is likely to be read by a specialist.
(Range of 20 — 39)




**Unacceptable quality
red**

If the image is lacking in quality, it is highly unlikely that a specialist will be able to read it.
(Range of 0 — 19)



About the Auto and Manual Exam modes

Auto Mode

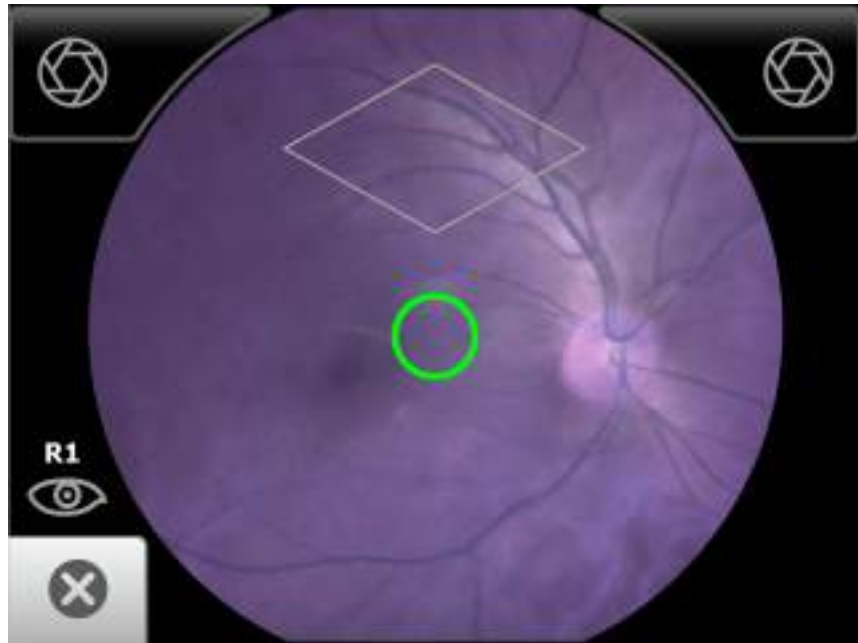


Auto Mode is the default image capture mode. Available automatic features include:

- Image capture
- Focus
- Flash brightness
- Sequential image capture order of the right and left eye
- Navigation to the image *Inspection screen*

In addition to the automatic features, manual capture, manual focus adjustment, and manual flash brightness adjustment are also available. (Manual focus adjustment and flash brightness are accessed from the *Capture settings screen*.)

Manual Mode

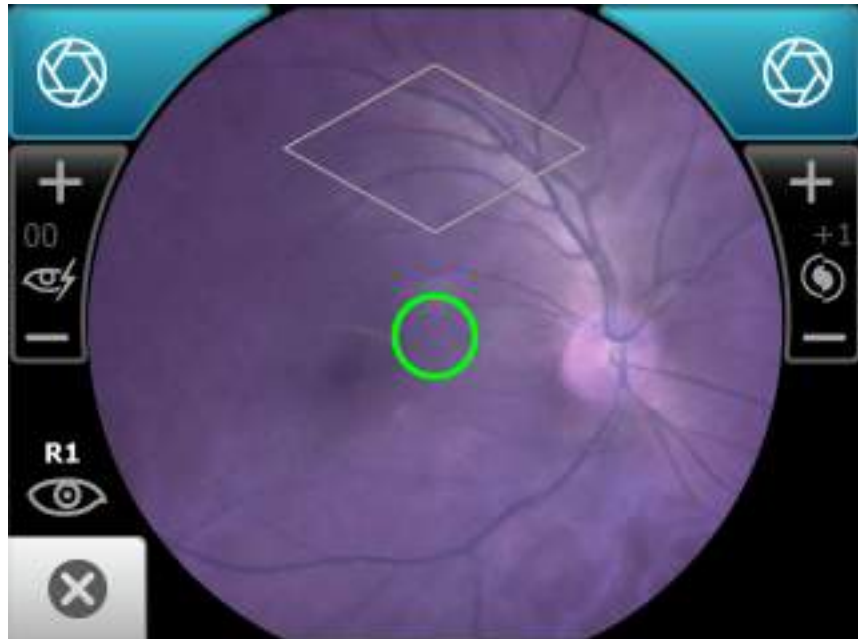

Manual mode contains advanced image capture options.

Note A manual touch of the capture button is required to obtain an image.


Available manual features include:

- Focus adjustment (values range from -1 to +1)
 - Flash brightness adjustment (values range from -2 to +2)
 - Sequential image capture order of the right and left eye
 - Image capture
-

Adaptive software in the Auto exam mode for software 6.10.02 and higher

Auto exam mode
first stage

Once the view is in-fundus, the adaptive software will enlarge the on-screen diamond at specific time intervals.

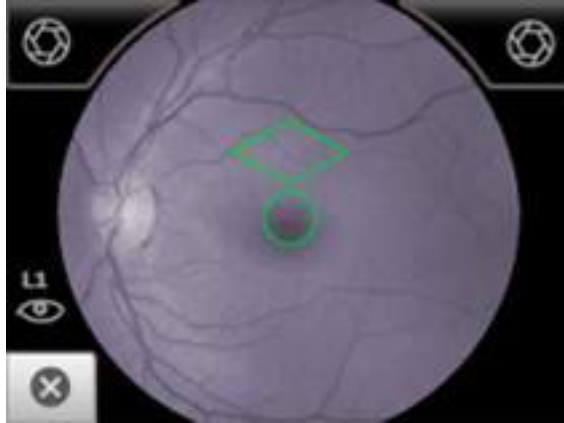

Note


The best image quality is obtained while the diamond is at the beginning smallest size.

Note

The diamond will reset to its smallest size each time the view moves out of the fundus.

Auto exam mode
second stage

As the image capture sequence continues, the diamond on the LCD touchscreen increases in size.

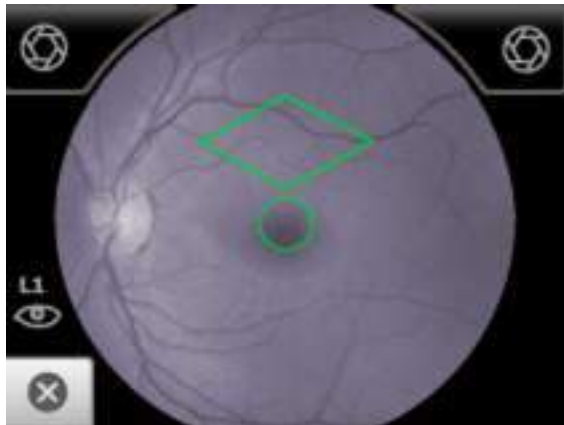

Perform an eye exam using the Auto exam mode

Auto exam mode is the default image capture mode.

Available automatic features include:

- Image capture
- Focus adjustment
- Flash brightness adjustment
- Sequential image capture order of the right and left eye
- Navigation to the image Inspection screen
 - In addition to automatic image capture, manual capture is also available.

Note To ensure that the patient's pupils sufficiently dilate to at least 3.5 mm diameter, adjust the room lighting to the lowest possible level. If necessary, have the patient sit in a dark room for 5 minutes to dilate their pupils.

Note While the RetinaVue 100 Imager can be used on patients with cataracts and other eye opacities, the use of the RetinaVue 100 Imager may result in a lower quality image due to the increased reflection of the flash off the patient's intraocular lens.

- Ensure that the SD card is installed into the RetinaVue 100 Imager.
- Ask your patient to remove their glasses, contacts can remain in place.
- Ensure that the patient sits on the edge of an exam table and stand in front of the patient to take the image.




Alternatively, ask the patient to sit in a chair and sit directly across from the patient with your legs together on the same side as the examined eye.



- Ask the patient to sit up straight and hold their head in a stationary position during the entire procedure.

- Ensure that the patient is sitting in a position that allows for a straight and level approach of the RetinaVue 100 Imager into the patient's eyes. A stable approach minimizes the loss of the view of the eye and repeated image capture attempts.
- Keep the RetinaVue 100 Imager lens the same height as the patient's examined eye.
- Approach the patient from directly in front of the examined eye.
- Make slight, micro adjustments with your guide hand to fine tune the position of the RetinaVue 100 Imager. Do not use your hand on the LCD side to make minor adjustments.

Note For additional information about techniques to capture the highest quality image, see the *Technique guide* and the *Training video* that came with the RetinaVue 100 Imager.


1. Press the power button .


Note To conserve power, the RetinaVue 100 Imager enters Sleep mode if it remains idle for approximately 2 minutes. Touch the LCD screen, press the power button, or press the Home button to wake up the RetinaVue 100 Imager from Sleep mode.

Note After 10 minutes of inactivity in Sleep mode, the RetinaVue 100 Imager powers down. Press the power button to activate the RetinaVue 100 Imager.

2. From the Home screen, touch **Exam**.


The *Preacquisition* screen appears.

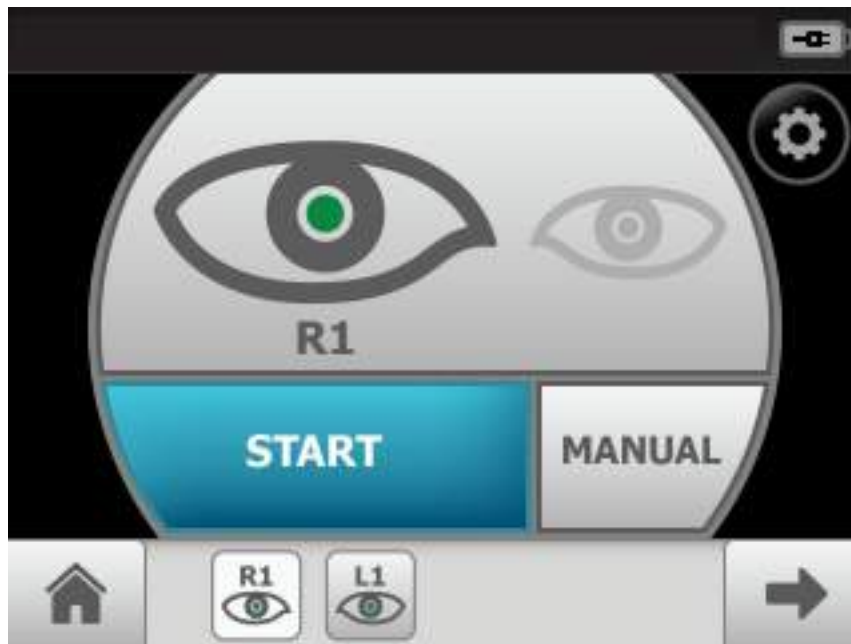

3. (Optional), If necessary, touch  (Settings) to adjust the illumination, or the focus. Touch the up or down arrows in the *Advanced image capture settings* and touch **OK**.

4. Touch **Start** to enter the *Image acquisition mode* and begin the first exam of the patient's right eye (R1).

The *Exam acquisition screen* appears.

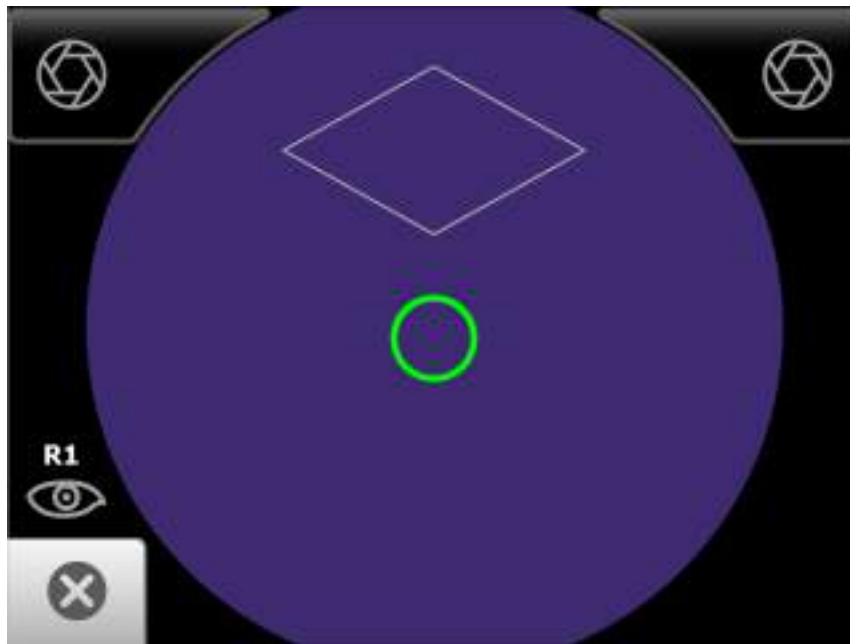

5. Hold the patient end of the RetinaVue 100 Imager two to three inches directly in front of the patient's examined eye. Continue forward to compress the eye cup against the examined eye.



WARNING Clean and disinfect the eye cup after each patient to avoid the risk of cross-contamination.




6. Direct the patient to focus on the green fixation lights inside the barrel of the RetinaVue 100 Imager.


Note Instruct the patient to cover, but not close, their unexamined eye. This will help the patient to focus on the green fixation lights.




7. Look for a white circular reflection at the top of the eye's image on the LCD touchscreen. This is your initial target.

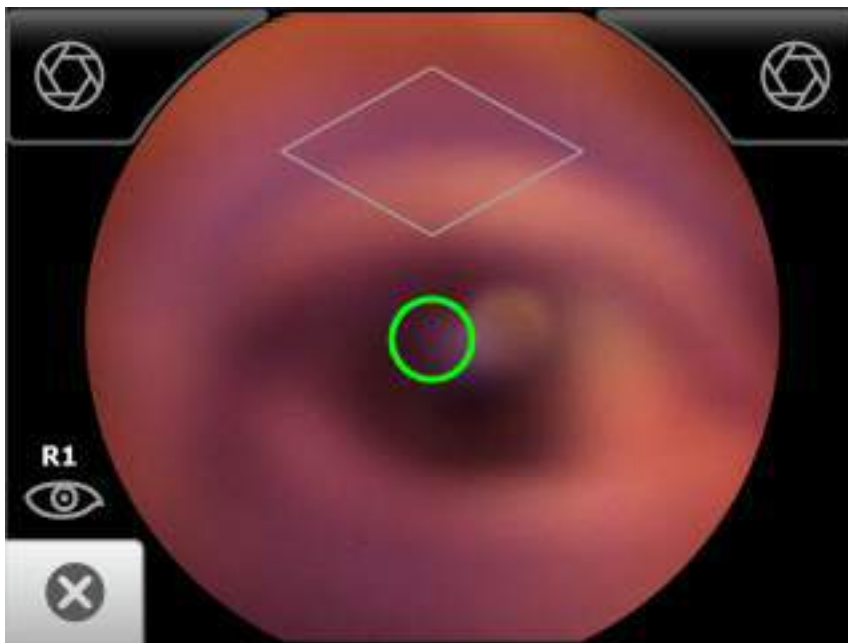


8. Guide the green circle on the RetinaVue 100 Imager LCD touchscreen towards the initial target by slowly moving the RetinaVue 100 Imager towards the eye.

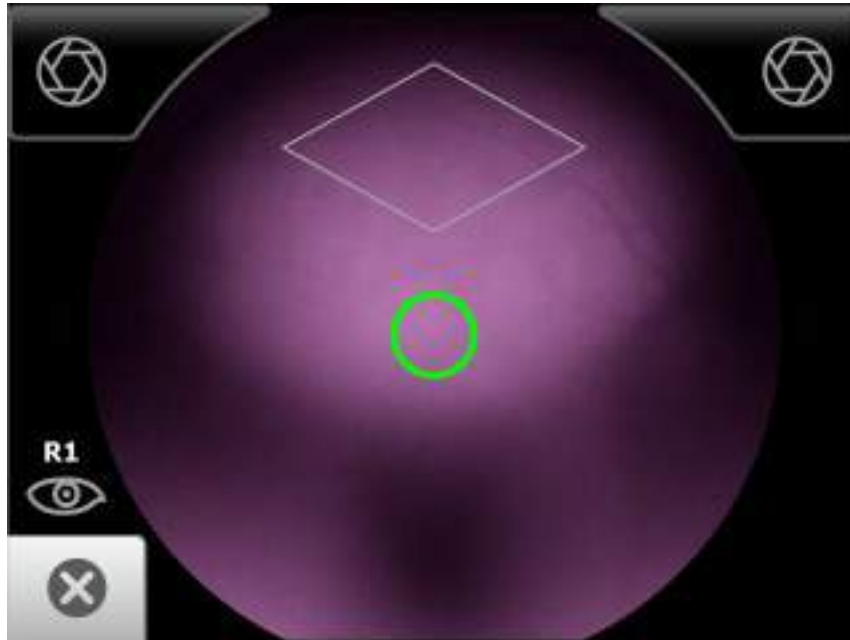

9. Once the reflection has filled the screen, push slightly forward until a “bubble” appears.
10. Keep the RetinaVue 100 Imager in line. Slightly adjust the RetinaVue 100 Imager away from any shadows that appear.
11. Move the patient end of the RetinaVue 100 Imager toward the bubble.

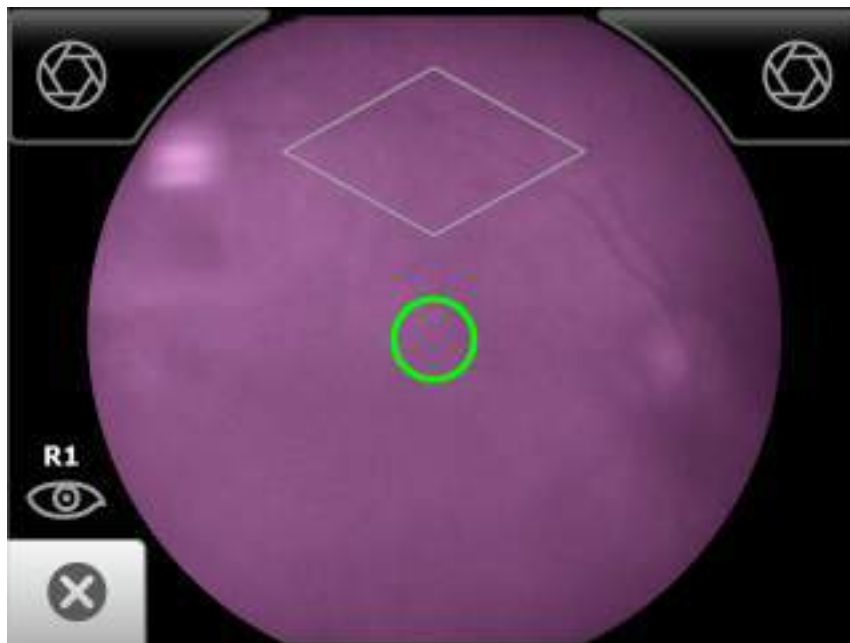


12. Align the bubble inside the diamond icon.

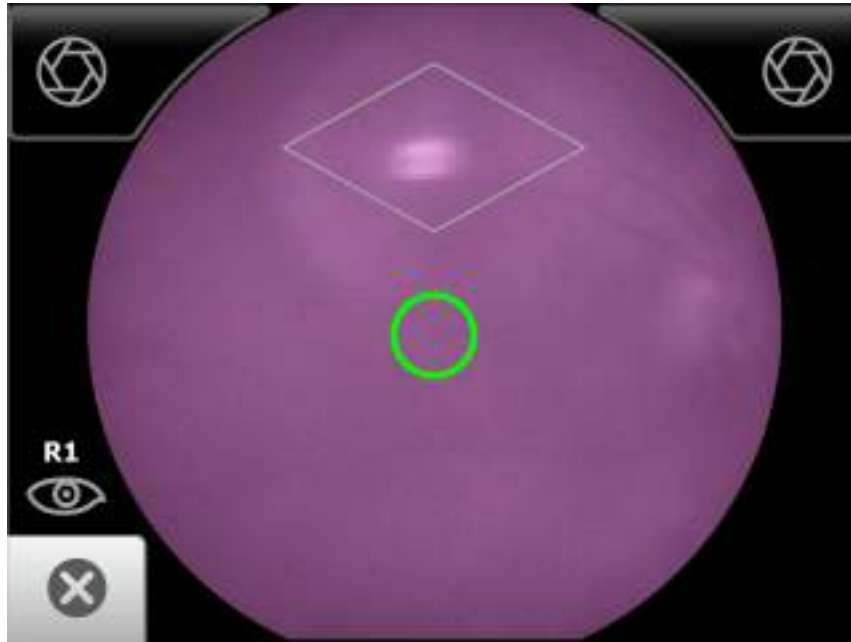






13. An audible tone alerts you that the image has been captured.


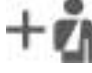
The acquired image appears on the LCD touchscreen for review in about five seconds.



14. Touch  (Accept) to begin the second exam of the patient's left eye (L1) and wait at least 1 minute before repeating the steps until the second image is automatically captured.
- (Optional) If the 2-image, Macula +Optic Disc-center exam setup is selected in Advanced settings, wait at least 1 minute before repeating the steps until all the images are automatically captured (R1, L1, R2, L2 and up to 4 additional images).



Note Touch  to review the current image. Images with a red X () denote low quality images. Touch  (navigates to the *Acquisition screen*) or touch  (navigates to the *Pre-acquisition screen*) to dispose of an image and then touch  (Add image) to begin a capture of another image.

- When the sequential image capture is complete and the images contain a green or yellow quality acceptance check mark, proceed to the *Exam summary screen*.


Note To add missing patient context information touch  and then touch  to create a new patient, or touch on an existing patient's name from within the *Patients screen* to select that patient. In the USB workflow, touch **Save** to save the exam images. In the wireless workflow, touch **Submit** to send the exam images to the RetinaVue Network Customer Portal.

Perform a manual image capture while in Auto exam mode

- From the *Pre-acquisition screen*, touch **Start**.

- Touch  or  (image capture) to manually capture an image.

An audible tone alerts you that the image has been captured. The acquired image appears on the LCD touchscreen for review in about five seconds.

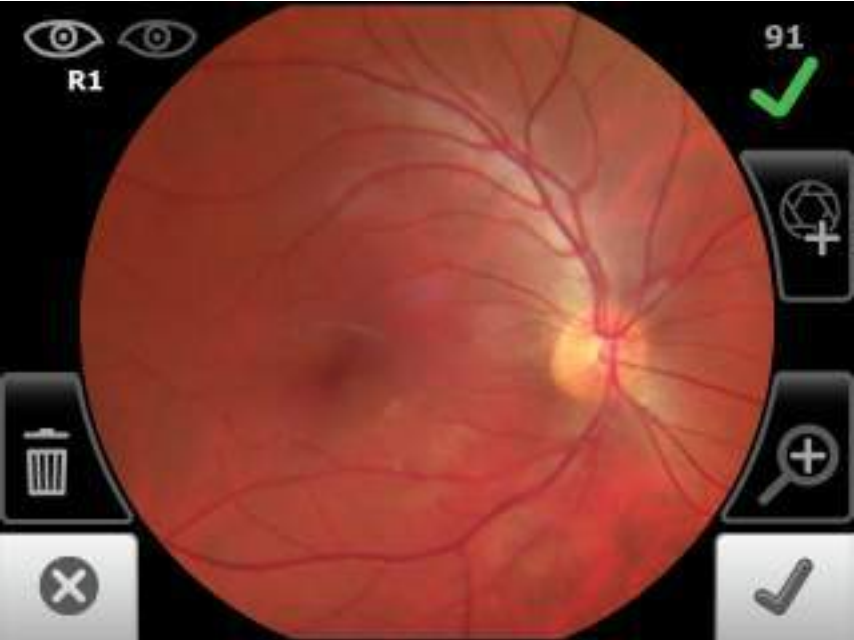






Skip an image capture while in Auto exam mode

1. From the Home screen, touch **Exam**.

The *Pre-acquisition screen* appears.




2. Touch  (Next) or touch , or , or  (Image Order) to skip to the next image in the sequence.
3. Touch **Start** to enter *Image acquisition mode* and begin the exam of the patient's eye using the next image in the sequence. The *Exam acquisition screen* appears.
4. If this is the last image in the sequence, touch  (Accept) on the *Image inspection screen* and navigate to the *Exam summary screen*. If an image is missing in the 1-image, Macula-centered Exam setup, touch **Override** to access the Save button (for USB) or the Submit button (for wireless) and to save or send the image. Or, touch **Override** to access the Save or Submit button and to save, or send, multiple images in the 2-image, Macula + Optic Disc-centered Exam setup.

Perform a manual image capture while in Manual exam mode


1. From the Home screen, touch **Exam**.



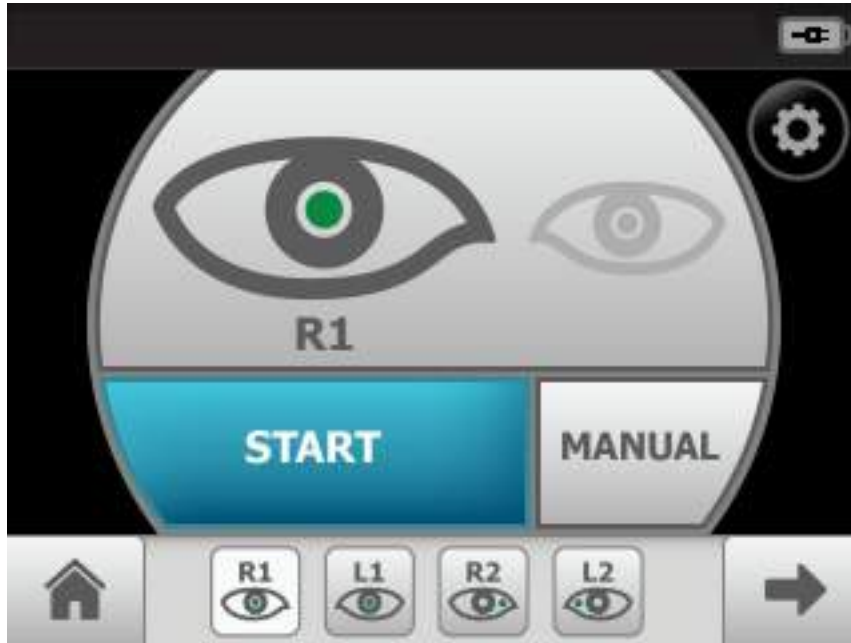

The *Pre-acquisition screen* appears.

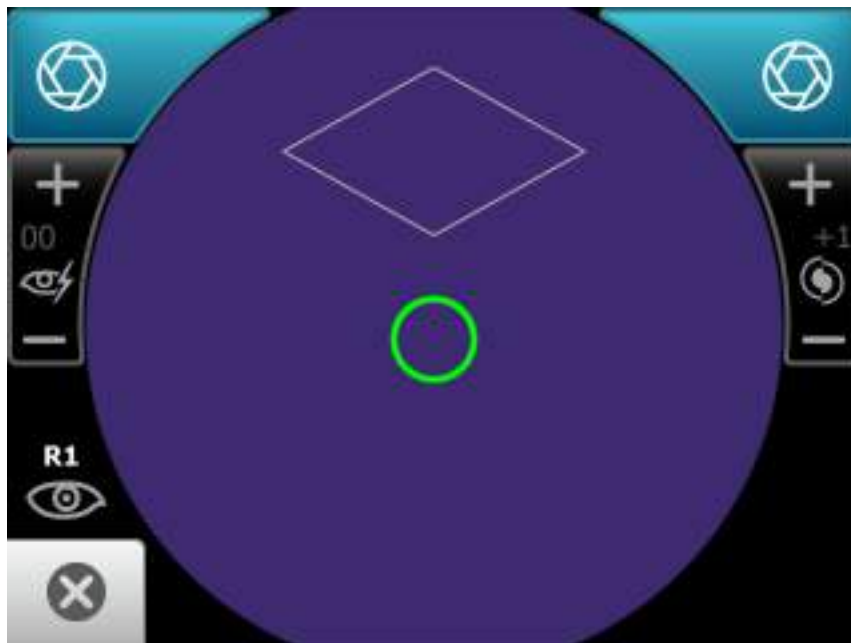


2. Touch **MANUAL** to enter *Image acquisition mode* and begin the exam of the patient's eye.

The *Exam acquisition screen* appears.



3. (Optional), If necessary, touch  to adjust illumination, or touch  to adjust the focus.
4. Hold the patient end of the RetinaVue 100 Imager two to three inches directly in front of the patient's examined eye. Continue forward to compress the eye cup against the examined eye.
5. Direct the patient to focus on the green fixation lights inside the barrel of the RetinaVue 100 Imager.
6. Look for a white circular reflection at the top of the eye's image on the LCD touchscreen. This is your initial target.
7. Guide the green circle on the RetinaVue 100 Imager LCD touchscreen towards the initial target by slowly moving the RetinaVue 100 Imager towards the eye.
8. Once the reflection has filled the screen, push slightly forward until a "bubble" appears.
9. Move the patient end of the RetinaVue 100 Imager toward the bubble.
10. Align the bubble inside the diamond icon.
11. Touch  or  (image capture) to manually capture an image.

An audible tone alerts you that the image has been captured. The acquired image appears on the LCD touchscreen for review in about five seconds.

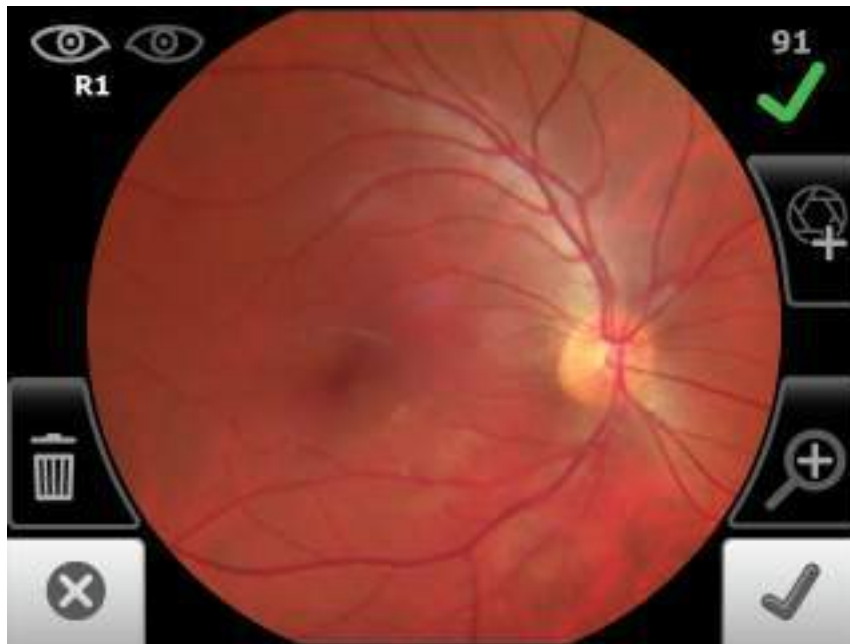

Troubleshooting




CAUTION A defective device should not be used. A device with parts which may be broken or missing, visibly worn, distorted, or contaminated should be replaced immediately.

To correct an error, follow these instructions or contact Welch Allyn Technical Support: www.welchallyn.com/service.

Troubleshooting image quality

Good quality image

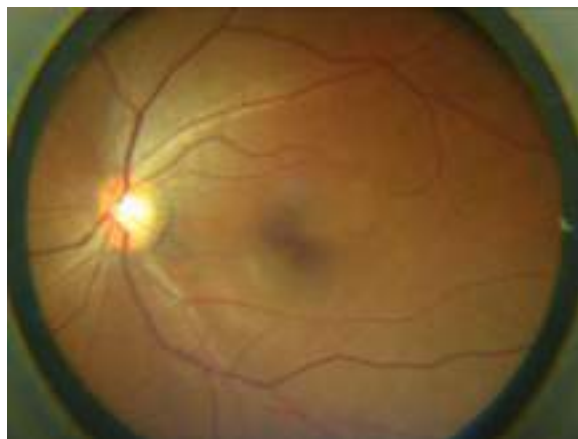

Problem

Partial blink




Solution

Ask the patient to keep his or her eye open.


Problem**Solution**

Dark semi-circle on top

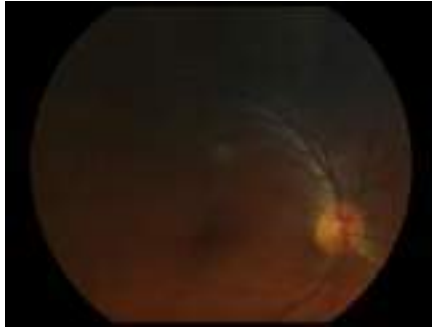

Move the patient side of the RetinaVue 100 Imager down and away from dark semi-circle.

Dirty lens

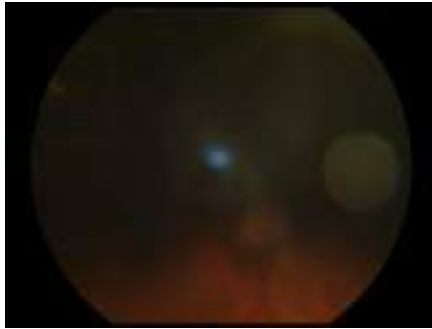

Clean the lens with an appropriate cleaning cloth.

Error messages

Error message/problem**Solution**

There is no SD card inserted. Please insert an SD card to continue. Insert SD card.

SD card access error

Check the write-protection switch on the SD card and ensure that it is in the forward position. Reinsert SD card to ensure the card is seated correctly.

The space left in the SD card is not enough. Please replace the SD card.

Replace the SD card with a formatted SD card.

Update file extraction failed. Power off device and restart.

Check the write-protection switch on the SD card and ensure that it is in the forward position.

Reinsert SD card to ensure the card is seated correctly.

If the SD card does not have sufficient space, remove all files (except the update image) from the SD card. Perform the update again.

Error message/problem	Solution
Failed saving files to the SD card.	Check the write-protection switch on the SD card and ensure that it is in the forward position. Reinsert or replace the SD card and click Retry.
Camera is docked. Remove camera from docking station to access camera functions.	Remove the camera from the docking station.
Place camera in charging dock with RetinaVue Network application running on the connected computer.	Connect your RetinaVue Imager charging dock to the PC running the RetinaVue Network application. Touch OK .
The battery is running low. Please recharge.	<p>When the battery charge level is below 25 percent, the low battery error message appears. Place the camera into the docking station to charge the battery.</p> <p>To reduce the likelihood of a low battery charging state, ensure that the USB cable is plugged into a "Always On" USB port on the computer.</p> <p>If the low battery message continues to display after the camera is docked in the docking station for several hours and while it is plugged into a USB port of your computer, plug the USB cable into an "Always On" USB port on the computer. (Only "Always On" USB ports continue to charge the camera when the computer enters sleep mode or power-saving mode.)</p>
The battery level is too low to acquire any image. Please recharge now.	When the battery charge level is below ten percent and the battery icon blinks, the low battery error message appears. Place the camera into the docking station to charge the battery.
Error connecting to the network.	Confirm that the network signal strength is sufficient. Re-enter the network name, username and password.
Error connecting to the network using the supplied credentials.	Re-enter the network name, username and password.

Problem**Solution**

Battery is not installed. Camera will not charge. The power indicator LED light on the camera blinks between green and amber when placed in the docking station.

Install the battery and check to ensure that the polarity on the battery matches the + and - polarity on the device.

Battery is installed with the polarity reversed. Camera will not charge. The power indicator LED light on the camera blinks between green and amber when placed in the docking station.

Re-install the battery and check to ensure that the polarity on the battery matches the + and - polarity on the device. Set the camera into the docking station and wait for the power indicator LED light to turn amber indicating that the device is charging.

The power indicator light does not turn on when the camera is set into the docking station. Charging is impacted if the pins are not properly touching.

Ensure that the contact pins of the camera are touching the contact pins of the docking station charger. If necessary, reseat the camera into the docking station.

The Windows® operating system may display the *Scan and fix* or *Continue without scanning* message once the camera is set in the cradle or when the USB cable is removed from the laptop or PC.

Select *Scan and fix* if this Windows message appears in situations when the camera is connected to a laptop or PC by the USB cable.



The Windows® operating system may display the AutoPlay notification popup dialogue box each time the camera is placed into the docking station when the USB cable is connected to the laptop or PC.

By default, the AutoPlay notification is displayed when a USB storage device is connected to the laptop or PC. To disable this notification, select **Control Panel > AutoPlay** from the Windows® Start menu. Uncheck the *Use AutoPlay for all media and devices* checkbox. Click **Save**.



Wireless connectivity

Problem

Cannot connect to the network.
Connection drops.

Solution

Ensure that the signal strength is sufficient. The signal strength appears in the top right corner of the device status bar.



Wireless signal strength

- Best (4 bars)
- Good (3 bars)
- Fair (2 bars)
- Weak (1 bar)
- No signal (no bars)
- No connection (X)

Ensure that Wi-Fi is connected.

Touch **Settings**.

Touch **Wi-Fi**.

Ensure that a network name appears. If no name appears, then select a network name again, enter a Username and Password and touch **OK**.

Cannot connect to the RetinaVue Network Server.
Message: *Error connecting to the network using the supplied credentials.*



Ensure that IP address is populated and shows a valid IP address other than 00:00:00:00:00:00.

Touch **Settings**.

Touch **Network Information**. The *Network Information* screen appears.

Wi-Fi MAC: A0:0A:BF1A:63:D7

AP MAC: C0:F8:DA:25:F3:6D

Wireless connectivity (Continued)

Problem	Solution
Clinic network does not appear in the Wi-Fi list.	<p data-bbox="964 352 1312 380">Level: 56% Noise 0 L.Q. 100</p> <p data-bbox="964 417 1211 445">IP Adr: 194.168.0.20</p> <p data-bbox="964 480 1240 508">Netmask: 266.266.266</p> <p data-bbox="964 543 1276 571">Broadcast: 194.168.0.255</p> <hr/> <p data-bbox="964 632 1393 693">Ensure that your network meets the supported credentials type.</p> <p data-bbox="964 728 1398 882">Check with your network administrator and ensure the network you are attempting to connect with is a Wi-Fi hot spot or has one of the following security types:</p> <p data-bbox="964 917 1357 978">EAP-PSK - pre-shared key using a passphrase</p> <p data-bbox="964 1014 1393 1075">EAP-PEAP - enterprise PEAP using a username and password</p> <p data-bbox="964 1110 1414 1199">There is no support for open/no security SSIDS or for certificate-based security types like EAP-TLS</p> <hr/> <p data-bbox="964 1293 1370 1320">Manually enter the network name.</p> <p data-bbox="964 1356 1154 1383">Touch Settings.</p> <p data-bbox="964 1419 1114 1446">Touch Wi-Fi.</p> <p data-bbox="964 1482 1406 1575">Scroll to the bottom of the list. Touch Other and enter the network name, followed by username and password.</p>

Wireless connectivity (Continued)

Problem

Solution



RetinaVue Imager camera does not connect to the network.

Ensure that the correct Username and Password are entered.

Touch **Settings**.

Touch **Wi-Fi**.

Select a network name from the list, enter your Username and Password, and check the "Remember this connection" checkbox. Every time you reboot the camera, it will connect to the network automatically.



Patient List is empty.

Follow all steps for connecting to RetinaVue Network Server.

Wireless connectivity (Continued)

Problem

Solution

Touch **Patients**.



Touch (refresh). Patients assigned to your clinic appear in the Patient list.



Follow all steps for connecting to RetinaVue Network Server.

Touch **Review**.

Any exams appearing in this list, have not been sent to the RetinaVue Network Server.

Exam was not sent to the RetinaVue Network Server.

Check your network connection. If your saved exam does not appear in the Review list, then it was successfully submitted to the RetinaVue Network Server. Check the "View Scheduled Patients" area on the RetinaVue Network Customer Portal. Select your clinic name from the drop-down menu. Find your patient and check the Status column. A "Submitted" status confirms that the exam was successfully received by the RetinaVue Network Server.



Wireless connectivity (Continued)**Problem****Solution**

The Save or Submit button is grayed out.

Ensure that you have a patient selected for the exam.

Ensure that you have the minimum number of exam images or have selected the override button to send fewer images than are required for the exam.

Manually entered Patients

For manually entered patients, ensure that the fields for Gender, Referring physician, and Diabetes type are filled in. If any of these fields are not filled in, leave the exam without saving the images and select from the *Home Screen/Patients*. Select a Gender using the radio buttons. Edit your patient's Referring Physician or Diabetes Type using the More button on the *Patient summary* screen. Begin the exam again by selecting OK and save new images for the exam.

Deleting pending exams from the *Restore factory defaults* menu

Problem**Solution**

In the event that you have an active Wi-Fi connection and unsent exams on the RetinaVue 100 Imager camera, the camera attempts to send exams before allowing the reset to factory defaults.



CAUTION Restoring the factory defaults will result in a loss of any images and patient data stored on the RetinaVue 100 Imager.

If you turn off Wi-Fi prior to attempting to reset factory defaults the camera will recognize that there is no connection available to send exams and will present a message about deleting data, however, the camera will still allow you to perform the factory default.

Touch **Settings**.

Touch **Wi-Fi On/Off**.

Touch **Advanced settings**.


Touch **Restore factory defaults**.

When the *Discard exams in Review and continue?* dialog appears, touch **Discard**.

When the *Reset all the settings to*

Problem	Solution
	<p><i>defaults?</i> dialog appears, touch Confirm.</p> <p>Touch OK.</p>

Reset the codes and data

Problem	Solution
<p>Locked out of the RetinaVue 100 Imager. Cannot recall the Device access code or the Advanced settings code.</p>	<div data-bbox="1019 653 1068 699" style="text-align: center;">  </div> <p data-bbox="1122 653 1421 768">CAUTION Resetting the codes will result in a loss of any images and patient data stored on the RetinaVue 100 Imager.</p> <p data-bbox="964 785 1421 974">When the <i>Please enter access code</i> screen appears, press and hold the Home button for approximately 10 seconds. Touch Confirm at the dialogue screen prompt <i>Reset all the settings to defaults?</i></p> <p data-bbox="964 1037 1421 1293">Note The reset process typically takes a few minutes to cycle through the complete reset. During the reset process, it is normal for the Confirm button to remain active before the reset confirmation dialogue screen appears. Do not remove the SD card or power off the device.</p> <p data-bbox="964 1314 1421 1404">Touch OK at the <i>All settings were restored to factory defaults</i> confirmation dialogue screen.</p>

Maintenance

Clean and disinfect



CAUTION Never submerge the RetinaVue 100 Imager in water or any other liquid.



CAUTION Never use abrasive cleaning agents, thinners or benzene for cleaning and never immerse the instrument in water or other cleaning liquids.



CAUTION Never insert a sharp object into the device or any open surface on the RetinaVue 100 Imager.



CAUTION Do not use unapproved cleaning or disinfection agents. Use of these agents may cause damage to components.



CAUTION When cleaning the device and the docking station charger, avoid using cloths or solutions that include quaternary ammonium compounds (ammonium chlorides) or glutaraldehyde-based disinfectants.



CAUTION Do not use bleach solutions of any kind when cleaning metal electrical contacts. They will damage the docking station charger and the device.

Note Clean on a routine basis according to your facility's protocols and standards or local regulations.

The RetinaVue 100 Imager has been tested to verify that it is compatible with the approved cleaning agents.

Cleaning the RetinaVue 100 Imager lens

Wipe the surface of the lens with a cleaning cloth or a lens cleaning tissue such as Thorlabs Inc.. Avoid touching the lens except when cleaning is required.

Cleaning the RetinaVue 100 Imager and eye cup

Power off the RetinaVue 100 Imager before cleaning it.

Use a pre-moistened wipe or soft cloth slightly moistened with an approved cleaner to gently clean the RetinaVue 100 Imager and eye cup.

Cleaning the docking station charger

Use a pre-moistened wipe or soft cloth slightly moistened with an approved cleaner to gently clean the docking station charger.

Disinfecting the RetinaVue 100 Imager, eye cup, and docking station charger

Note Disinfect according to your facility's protocols and standards or local regulations.

Use a pre-moistened wipe or soft cloth slightly moistened with an approved disinfectant to disinfect the RetinaVue 100 Imager, eye cup, and docking station charger. Follow the manufacturers recommendations.

Approved cleaning agents

Use only the following approved cleaning agents on the RetinaVue 100 Imager, eye cup, and docking station charger:

- 70 percent isopropyl alcohol
- CaviWipes
- Sani-Cloth® Plus
- Super Sani-Cloth®

The RetinaVue 100 Imager has been tested to verify these instructions as being capable of preparing the RetinaVue 100 Imager for re-use. You must ensure that cleaning and disinfection as actually performed by your personnel, with your equipment and materials, achieves the desired result. This may require validation and routine monitoring of your actual process.

Specifications

Specifications

View Angle	45 Degree (Typical)
Diopter	-20 ~ +20D (Typical)
Dimensions	L 8.8 inches (224.2 mm) × W 4.2 inches (106.7 mm) × H 4.1 inches (105.2 mm)
Weight	565 Grams ±5%
Search Fundus Lighting	Infrared LED
Camera flash light	Natural White Light Emitting Diode (LED)
Captures	Auto capture and Manual capture
Camera Resolution	2048 x 1536 pixels
LCD Monitor	3.5 inch Full Color TFT-LCD
Image Format	Encrypted JPEG/DICOM
Interface	Micro USB
File Transfer	Docking station, micro USB cable to PC running RetinaVue Network (requires additional software)
File Storage	SD Card, default 8GB. Supports 8GB to 32GB by FAT32 format. Supports 1,200 exams Supports up to 200 patients
Power Source	Rechargeable Lithium Battery 3.6V / 3350mAh
External Power	Source: 100~240 VAC, 50/60 Hz
Power Adapter Spec.	Input Spec. : 100~240 VAC, 50/60 Hz, 0.2A - 0.1A; Output Spec. : 5V DC, 2A
Charger Station Input Spec.	5V DC, 1.2A

Operating Time	3 hours at 5 minutes working mode. The device enters Sleep mode if idle for 2 minutes. Once the device is idle for 2 minutes it enters Sleep mode and if it remains idle in Sleep mode for another 10 minutes the device powers down.
----------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

Charging Time	5 hours (through the docking station using a Welch Allyn approved battery)
---------------	----------------------------------------------------------------------------

All specifications are subject to change.

Operating environment

Criterion	Environmental conditions
Temperature	+ 10 °C to + 35 °C
Relative humidity	30 % to 90 %
Atmospheric pressure	800 hPa to 1060 hPa

Storage environment

Criterion	Environmental conditions
Temperature	- 10 °C to + 55 °C
Relative humidity	10 % to 95 %
Atmospheric pressure	700 hPa to 1060 hPa

Transportation environment

Criterion	Environmental conditions
Temperature	- 40 °C to + 70 °C
Relative humidity	10 % to 95 %
Atmospheric pressure	500 hPa to 1060 hPa

Device lifetime

The lifetime of the device is five (5) years. The manufacture date of the device can be found on the device labels.

Welch Allyn will service RetinaVue 100 Imager devices that are within their lifetime. Firmware updates and support may require an annual subscription service after the initial one (1) year warranty period. The expected battery life is at least one (1) year. If the RetinaVue 100 Imager device fails to hold a charge, a new battery can be ordered.

Device radio

The RetinaVue 100 Imager's radio operates on 802.11 networks.

WLAN standard	IEEE802.11b/g/n, Wi-Fi compliant
Frequency Range	2.412 to 2.484 GHz
Data Rate	802.11b: 11, 5.5, 2, 1 Mbps DSSS 802.11g: 54, 48, 36, 24, 18, 12, 9, 6 Mbps OFDM 802.11n: HT20 MCS0~7, HT40 MCS0~7
Modulation	802.11b: CCK, DQPSK, DBPSK 802.11g: 64 QAM, 16 QAM, QPSK, BPSK 802.11n: BPSK, QPSK, 16-QAM, 64-QAM
Host interface	USB 2.0
Encryption	WPA2, WPA2-PSK, WPA2-Enterprise (EAP-PEAP), IEEE 802.11i
Channel availability	2.4GHz: (Ch. 1-11) – United States (Ch. 1-13) – Europe (Ch. 1-14) – Japan
Antenna	PCB-37x6-15U.FL with I-PEX cable, Max Gain: 1.3 dBi
Agency approvals	US: FCC 15.247:2014, KDB 558074 V3, Part 18, ANSI C63.10.2013 Europe: Annex 3.2 RED in 2.4GHz bands, EN 300 328 V. 2.1.1, EN 301 489-17 V3.2.0, EN 301 489-1 V2.2.0 and EN 62479

Australia and New Zealand: Australian Communications and Media Authority (ACMA) Radio Compliance Mark (RCM)



New Zealand maintains a Mutual Recognition Agreement (MRA) with Australia.

To ensure compliance with local regulations, be sure the correct country in which the access point is installed is selected. This product can be used with the following restriction(s):

Norway - Does not apply for the geographical area within a radius of 20 km from the center of Ny-Ålesund.

France - Outdoor use is limited to 10 mW EIRP within the band 2454 to 2483.5 MHz.

Note Effective Isotropic Radiated Power (EIRP).

General compliance and standards

The RetinaVue 100 Imager complies with the following standards:

IEC 60601-1:2005/A1:2012 (EN 60601-1:2006/A1:2013)

IEC 60601-1-2 : 2014 (EN 60601-1-2 : 2015)

ISO 15004-1 : 2006

ISO 15004-2 : 2007

ISO 10940 : 2009

ISO 10993-1 : 2009



**Directive 2002/96/EC-WEEE:
Disposal of noncontaminated electrical and electronic equipment**

This product and its accessories must be disposed of according to local laws and regulations. Do not dispose of this product as unsorted municipal waste. Prepare this product for reuse or separate collection as specified by Directive 2002/96/EC of the European Parliament and the Council of the European Union on Waste Electrical and Electronic Equipment (WEEE). If this product is contaminated, this directive does not apply.

For more specific disposal or compliance information, see www.welchallyn.com/weee, or contact Welch Allyn Customer Service at +44 207 365 6780.

General radio compliance

The wireless features of this device must be used in strict accordance with the manufacturer's instructions as described in the user documentation that comes with the product.

Federal Communications Commission (FCC)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- This device may not cause harmful interference.
- This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy. If not installed and used in accordance with the instructions, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try and correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the distance between the equipment and the receiver
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected
- Consult the dealer or an experienced radio/TV technician for help

The user may find the following booklet prepared by the Federal Communications Commission helpful:

The Interference Handbook

This booklet is available from the U.S. Government Printing Office, Washington, D.C. 20402. Stock No. 004-000-0034504.

Welch Allyn is not responsible for any radio or television interference caused by unauthorized modification of the devices included with this Welch Allyn product, or the substitution or attachment of connecting cables and equipment other than specified by Welch Allyn.

The correction of interference caused by such unauthorized modification, substitution, or attachment will be the responsibility of the user.

The device is designed and tested to meet the applicable limits for radio frequency (RF) exposure established by the Federal Communications Commission (U.S.A.) Specific Absorption Rate; (SAR) refers to the rate at which the body absorbs RF energy.

The maximum allowable SAR limit is 4.0 W/kg, averaged over 10 gram of tissue for the device. The following SAR values were measured during product certification, 1.060 W/kg with a (0mm gap) a value well below the maximum allowable limit of 4.0 W/kg.

Industry Canada (IC) emissions

This device complies with Industry Canada's licence-exempt RSSs.

Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of this device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence.

L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage; (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

This Class B digital apparatus complies with Canadian ICES-003, CAN ICES-3(B)/NMB-3(B).

Cet appareil numérique de classe B est conforme à la norme ICES-003 du Canada, CAN ICES-3 (B) / NMB-3 (B).

RF Radiation Hazard Warning

Using higher gain antennas and types of antennas not certified for use with this product is not allowed. The device shall not be co-located with another transmitter.

Cet avertissement de sécurité est conforme aux limites d'exposition définies par la norme CNR-102 at relative aux fréquences radio.

This radio transmitter (Contains IC ID: 21375-WAWEC100) has been approved by Industry Canada to operate with the antenna types listed in table above with the maximum permissible gain and required antenna impedance for each antenna type indicated. Antenna types not included in this list, having a gain greater than the maximum gain indicated for that type, are strictly prohibited for use with this device.

Le présent émetteur radio (Contains IC ID: 21375-WAWEC100) a été approuvé par Industrie Canada pour fonctionner avec les types d'antenne énumérés ci-dessous et ayant un gain admissible maximal et l'impédance requise pour chaque type d'antenne. Les types d'antenne non inclus dans cette liste, ou dont le gain est supérieur au gain maximal indiqué, sont strictement interdits pour l'exploitation de l'émetteur.

European Union

Czech	Medimaging Integrated Solution Inc. (MiiS) tímto prohlašuje, že tento vysílač s nízkým výkonem je v souladu se základními požadavky a dalšími příslušnými ustanoveními směrnice 2014/53/ ES
Danish	Medimaging Integrated Solution Inc. (MiiS) erklærer herved, at denne lavt strøm transmitter er i overensstemmelse med de væsentlige krav og andre relevante bestemmelser i direktiv 2014/53/EF
Dutch	Hierbij verklaart Medimaging Integrated Solution Inc. (MiiS) dat deze low power transmitter voldoet aan de essentiële eisen en andere relevante bepalingen van Richtlijn 2014/53/EG
English	Hereby, Medimaging Integrated Solution Inc. (MiiS) declares that this low power transmitter is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EC.
Estonian	Käesolevaga Medimaging Integrated Solution Inc. (MiiS) teatab, et see väikese võimsusega saatja on vastavus olulistele nõuetele ja teistele asjakohastele sätetele direktiivi 2014/53/ EÜ
Finnish	Näin ollen Medimaging Integrated Solution Inc. (MiiS) vakuuttaa, että tämä pienitehoinen lähetin on direktiivin 2014/53/EY olennaisten vaatimusten ja muiden asiaa koskevien säännösten mukainen
French	Par conséquent, Medimaging Integrated Solution Inc. (MiiS) déclare que cet émetteur de faible puissance est conforme aux exigences essentielles et autres dispositions pertinentes de la directive 2014/53/CE
German	Hiermit erklärt Medimaging Integrated Solution Inc. (MiiS), dass dieser Niedrigleistungssender den grundlegenden Anforderungen und anderen relevanten Bestimmungen der Richtlinie 2014/53/EG entspricht
Greek	Με αυτό τον τρόπο, η Medimaging Integrated Solution Inc. (MiiS) δηλώνει ότι αυτός ο πομπός χαμηλής ισχύος είναι σύμφωνος με τις βασικές απαιτήσεις και άλλες σχετικές διατάξεις της οδηγίας 2014/53/EK

Hungarian	A Medimaging Integrated Solution Inc. (MiiS) ezúton kijelenti, hogy ez az alacsony teljesítményű adó megfelel a 2014/53/EK irányelv alapvető követelményeinek és egyéb vonatkozó rendelkezéseinek
Italian	In questo modo, Medimaging Integrated Solution Inc. (MiiS) dichiara che questo trasmettitore di bassa potenza è conforme ai requisiti essenziali e alle altre pertinenti disposizioni della direttiva 2014/53/CE.
Latvian	Ar šo, Medimaging Integrated Solution Inc. (MiiS) paziņo, ka šis mazjaudas raidītājs atbilst būtiskajām prasībām un citiem attiecīgiem noteikumiem Direktīvā 2014/53/EK
Lithuanian	Šiuo dokumentu Medimaging Integrated Solution Inc. (MiiS) deklaruoja, kad šis mažas energijos siųstuvas atitinka esminius reikalavimus ir kitas susijusias nuostatas Direktyva 2014/53/EB
Malti	Hawnhekk, Medimaging Integrated Solution Inc. (MiiS) jiddikjara li din trasmettitur enerġija baxxa hija konformi mar-rekwiziti essenzjali u dispozizzjonijiet rilevanti oħra tad-Direttiva 2014/53/KE
Polish	Niniejszym Medimaging Integrated Solution Inc. (MiiS) oświadcza, że ten nadajnik o małej mocy jest zgodny z zasadniczymi wymaganiami i innymi istotnymi przepisami dyrektywy 2014/53/ WE
Portuguese	Por este motivo, a Medimaging Integrated Solution Inc. (MiiS) declara que este transmissor de baixa potência está em conformidade com os requisitos essenciais e outras disposições relevantes da Diretiva 2014/53/CE
Slovak	Medimaging Integrated Solution Inc. (MiiS) týmto vyhlasuje, že tento vysielač s nízkym výkonom je v súlade so základnými požiadavkami a ďalšími príslušnými ustanoveniami smernice 2014/53/ES
Slovene	S tem, Medimaging Integrated Solution Inc. (MiiS) izjavlja, da je ta nizka moč oddajnik v skladu z bistvenimi zahtevami in drugimi ustreznimi določbami Direktive 2014/53/ES
Spanish	Por este medio, Medimaging Integrated Solution Inc. (MiiS) declara que este transmisor de baja potencia cumple con los requisitos esenciales y otras disposiciones pertinentes de la Directiva 2014/53/EC
Swedish	Medimaging Integrated Solution Inc. (MiiS) förklarar härigenom att denna låg effekt sändare överensstämmer med de väsentliga kraven och andra relevanta bestämmelser i direktiv 2014/53/EG

EMC compliance

Special precautions concerning electromagnetic compatibility (EMC) must be taken for all medical electrical equipment. This device complies with IEC 60601-1-2 : 2014 (EN 60601-1-2 : 2015).

- All medical electrical equipment must be installed and put into service in accordance with the EMC information provided in this *Directions for use*.
- Portable and mobile RF communications equipment can affect the behavior of medical electrical equipment.


The 901103 Imager complies with all applicable and required standards for electromagnetic interference.

- It does not normally affect nearby equipment and devices.
- It is not normally affected by nearby equipment and devices.
- However, it is good practice to avoid using the 901103 Imager in extremely close proximity to other equipment.
- It is not safe to operate the 901103 Imager in the presence of high-frequency surgical equipment.

Emissions and immunity information

Electromagnetic emissions

The 901103 is intended for use in the electromagnetic environment specified below. The customer or user of the 901103 should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The 901103 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The 901103 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	 WARNING This equipment/system is intended for use by healthcare professionals only. This equipment/ system may cause radio interference or may disrupt the operation of nearby equipment ^a . It may be necessary to take mitigation measures, such as re-orienting or relocating the 901103 or shielding the location.

^a The 901103 contains a 5-GHz orthogonal frequency-division multiplexing transmitter or a 2.4-GHz frequency hopping spread-spectrum transmitter for the purpose of wireless communication. The radio is operated according to the requirements of various agencies, including FCC 47 CFR 15.247 and R&TTE Directive (1995/5/EC). The transmitter is excluded from the EMC requirements of 60601-1-2, but should be considered when addressing possible interference issues between this and other devices.

Electromagnetic immunity

The 901103 is intended for use in the electromagnetic environment specified below. The customer or the user of the 901103 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	± 1 kV differential mode ± 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on	100% dip in 0.5 cycle 100% dip in 1 cycles 30% dip in 25 cycles	100% dip in 0.5 cycle 100% dip in 1 cycles 30% dip in 25 cycles	Mains power quality should be that of a typical commercial or hospital environment. If the user of the 901103 requires continued operation during power mains interruptions, it is

Electromagnetic immunity

power supply input lines IEC 61000-4-11	100% interruption in 250 cycles	100% interruption in 250 cycles	recommended that the 901103 be powered from an uninterruptible power supply or a battery.
--------------------------------------------	---------------------------------	---------------------------------	-------------------------------------------------------------------------------------------

Electromagnetic immunity

The 901103 is intended for use in the electromagnetic environment specified below. The customer or the user of the 901103 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
			Portable and mobile RF communications equipment should be used no closer to any part of the 901103, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Recommended separation distance			
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	$d = (1.17) \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 1 GHz	3 V/m	$d = (1.17) \sqrt{P}$ 80 to 800 MHz

$$d = (2.33) \sqrt{P} \text{ 800 MHz to 2,7 GHz}$$

where P is the maximum output power rating of the transmitter in watts (W) and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol:



Note1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

^aField strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the 901103 is used exceeds the applicable RF compliance level above, the 901103 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the 901103.

^bOver the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the 901103

The 901103 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or user of the 901103 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the 901103 as recommended below, according to the maximum output power of the communications equipment.

Separation distance according to frequency of transmitter (m)

Rated max. output power of transmitter (W)	150 kHz to 80 MHz $d = (1.17) \sqrt{P}$	80 MHz to 800 MHz $d = (1.17) \sqrt{P}$	800 MHz to 2.7 GHz $d = (2.23) \sqrt{P}$
0.01	0.11667	0.11667	0.23333
0.1	0.36894	0.36894	0.73785
1	1.1667	1.1667	2.3333
10	3.6894	3.6894	7.3785
100	11.667	11.667	23.3333

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.


Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

Appendix

RetinaVue 100 Imager Home screen

The RetinaVue 100 Imager *Home screen* includes the following areas:

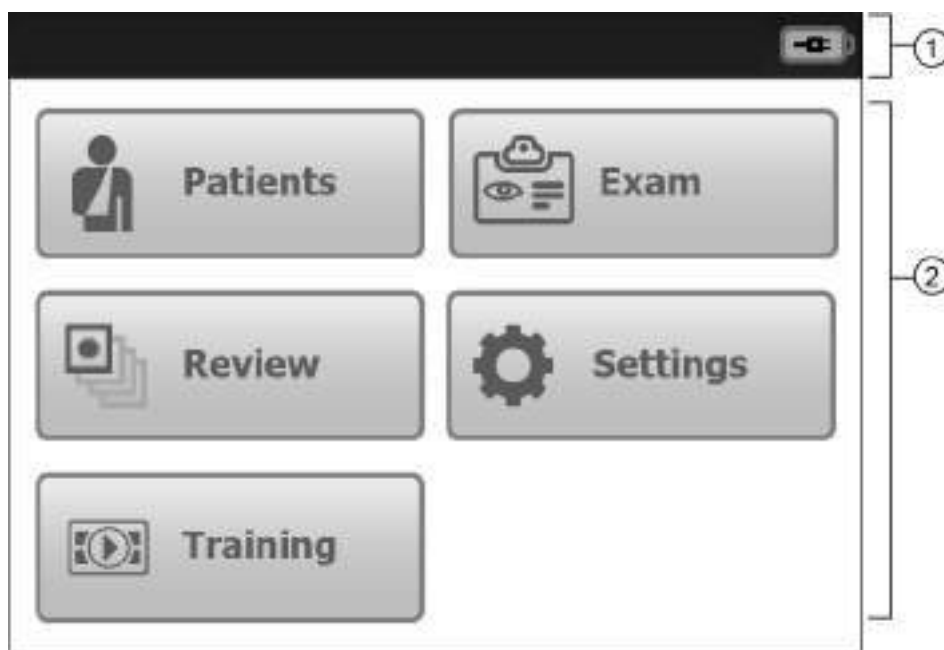

Item	Area
1	Device status
2	Content




Device Status Area

The Device Status Area, located at the top of the *Home screen*, displays the Battery condition. Estimated battery capacity is displayed in a charged or discharged status format.

In the *Exam Summary* and *Review screens*, the Device Status Area displays:



- Patient Icon and Patient name. The format of the Patient name is last name, first name, and middle initial.
- Connection status. The icons indicate which connection type, if any, is currently active.

Icon	Connection type
	Wireless signal strength
	Expected connection to the the RetinaVue Network Server is not occurring.
	USB

- Battery condition. Estimated battery capacity is displayed in a charged or discharged status format.

Content Area

The Content Area includes four workflow selection buttons and a Training button. Touching one of the five buttons navigates to the next screen associated with that workflow or training.

- **Patients:** Presents a list of patients created on the RetinaVue 100 Imager or transferred from the RetinaVue™ Network (requires additional software).
- **Exam:** The Exam workflow provides a guided approach to image capture. To maintain consistent and complete exam submissions, a prescribed image acquisition sequence is used to complete an exam. The Exam workflow also presents controls, status, and previews of images. (The Exam Workflow requires additional RetinaVue™ Network software to submit images to an overread service.)
- **Review:** Presents images and exams from a list of patients contained on the RetinaVue 100 Imager.
- **Settings.** Provides controls to view or change the following:
 - Wi-Fi On/Off
 - Wi-Fi
 - Clinic /States
 - Volume
 - Time
 - Date
 - Advanced settings:
 - Device access code (USB only)
 - Advanced settings code
 - Format SD
 - Restore factory defaults
 - Exam setup
 - Export log

- About
 - Software version
- Regulatory Information
- Network Information
- **Training:** Provides demonstrations and interactive training exercises on the RetinaVue 100 Imager. Short video clips and interactive training help you practice the three key steps necessary to master the RetinaVue 100 Imager.

Navigation controls in the Training mode


From the *Home screen* touch **Training**. The RetinaVue 100 Imager *Interactive Training Video screen* appears.


Note Navigation controls change with each mode.


Touch anywhere within the RetinaVue 100 Imager touchscreen to display the training navigation controls. The navigation controls appear within the middle of the RetinaVue 100 Imager touchscreen.




The navigation controls for the Interactive Training video include:


Step back to beginning: Touch  to navigate to the first scene within the training section.

Step back to previous scene: Touch  to navigate to the previous scene within the training section.

Home:  Displays the RetinaVue 100 Imager workflow selection buttons and navigates to the *Home screen*.

Play: Touch  to resume the training.

Skip forward to next scene: Touch  to navigate to the next scene within the training section.

Skip forward to last: Touch  to navigate to the last scene within the training section.

Note In addition to the control buttons, touch anywhere within the RetinaVue 100 Imager touchscreen to pause the training or to display the training navigation controls.

Navigation controls in the Exam mode

The Navigation area is located at the bottom of the RetinaVue 100 Imager screen. Navigation controls change with the mode. The navigation controls for Pre-, post-, and acquisition in Exam mode include:

Pre-acquisition

1-image


Macula-centered



2-image

Macula + Optic Disc-centered



Home:  Displays the RetinaVue 100 Imager workflow selection buttons and navigates to the *Home screen*.







Next:  Navigates to the next image acquisition within a prescribed image acquisition sequence or to the next step in a sequence of steps.

Image acquisition navigation: Navigates to the image acquisition position within the prescribed sequence. Provides a visual representation of the prescribed sequential

acquisition, including the right eye ( ) and the left eye ( ). Provides a status of acquired images for acceptance or dismissal.

Acquisition

Exit:  Navigates to the *Pre-acquisition screen*.

Post-acquisition

Discard image:  Navigates to the *Pre-acquisition screen*. Or  Navigates to the *Acquisition screen*.




Accept image:  Navigates to the *Exam summary screen* when the last exam in the sequence has been accepted. Also navigates to *Pre-acquisition screen* if there are more images to be captured.


Image inspection:  Navigates to the *Image inspection screen*.

Add image:  Navigates to the *Acquisition screen*

Navigation controls in the Settings mode

The navigation controls for the Settings mode include the back, home, and advance to next buttons. The primary navigation area is located at the bottom of the RetinaVue 100 Imager screen. While on the *Settings* screen, the > (Advance to next screen) button navigates to the next dialogue. The button is located on the same line as the setting.


Back:  Navigates to a previous dialogue or step.


Home:  Displays the RetinaVue 100 Imager workflow selection buttons and navigates to the *Home screen*.

>: (Advance to next screen) Navigates to a next dialogue while on the *Settings* screen.


Navigation controls in the Review mode

The navigation controls for the Review mode include the back, home, advance to next screen, and image quality buttons. The primary navigation area is located at the bottom of the RetinaVue 100 Imager screen. In addition, the > (Advance to next screen) button navigates to an image. The button is located on the same line as the patient's name on the *Review* screen. The image quality button also navigates to an image. The button is located below the patient's name on the *Exam summary* screen.

Back:  Navigates to a previous dialogue or step.

Home:  Displays the RetinaVue 100 Imager workflow selection buttons and navigates to the *Home screen*.

>: (Advance to next screen) Navigates to an image while on the *Review* screen.

Good quality green:  From the *Exam summary* screen, touch the good quality icon to navigate to the image on the *Post-acquisition* screen.

Acceptable quality yellow:  From the *Exam summary* screen, touch the acceptable quality icon to navigate to the image on the *Post-acquisition* screen.


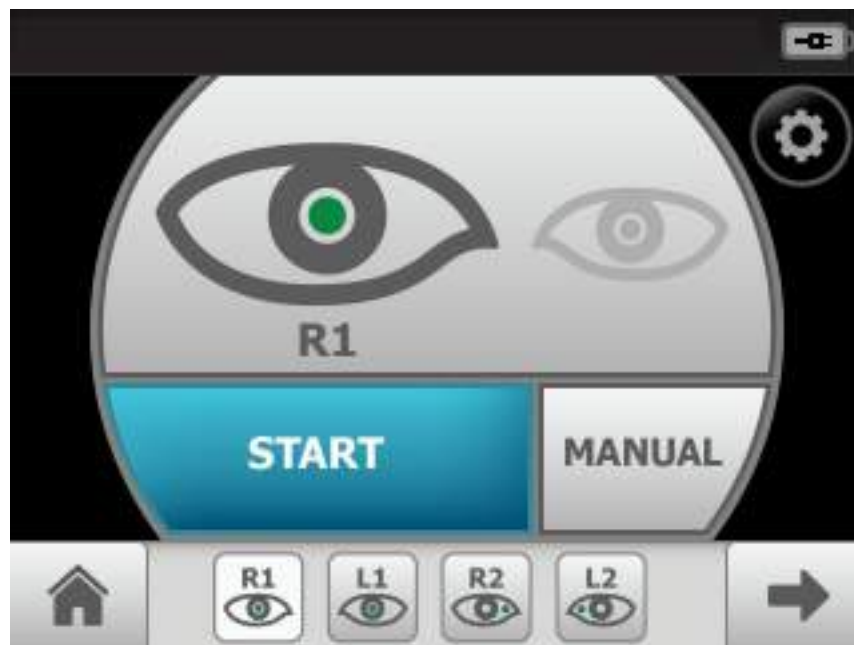




Unacceptable quality red:  From the *Exam summary* screen, touch the unacceptable quality icon to navigate to the image on the *Post-acquisition* screen.




Image Pre-acquisition controls (Auto mode)


From the *Home screen* touch **Exam**. The *Exam Pre-acquisition screen* appears.







No.	Feature	No.	Feature
1.	START enters automatic image acquisition mode.	4.	 Navigates to the <i>Home screen</i> .
2.	 (Image Acquisition Settings) navigates to the <i>Capture settings</i> screen. Adjust the flash brightness (-2 to +2) or the focus settings (-1 to +1).	5.	MANUAL switches to advanced image capture options.
3.	 (Image Order) skips to the next image in the sequence.	6.	 (Next) skips to the next image acquisition within a prescribed image acquisition sequence or navigates to the <i>Exam summary screen</i> if the last image has been acquired.

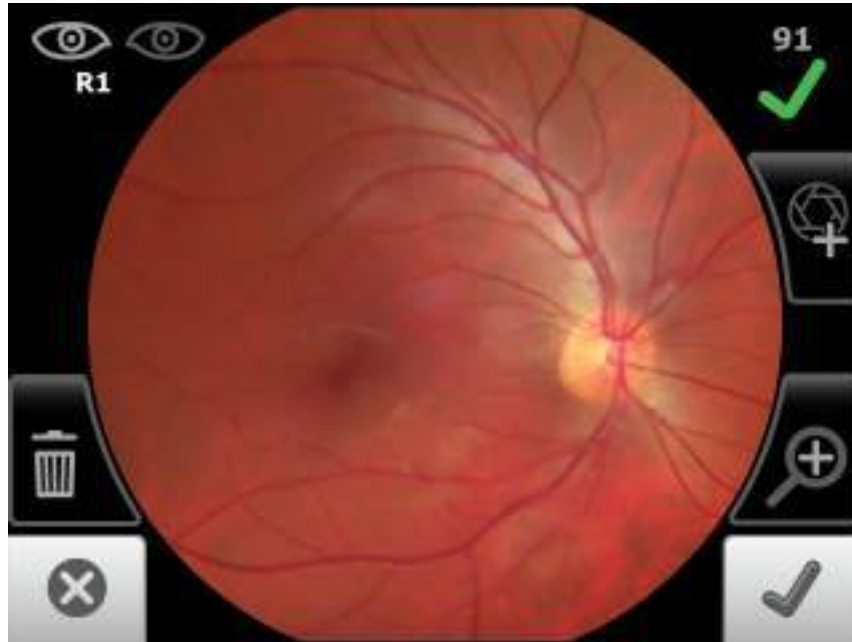
Controls (Post-acquisition)

On the *Post-acquisition screen*, images with a green check () are of good quality (scores of 40 +), images with a yellow check () are of acceptable quality (scores of 20 — 39), and those with a red x () are of unacceptable quality (scores of 0 — 19).





Use the RetinaVue 100 Imager touchscreen inspection control  to navigate to the *Inspection screen*, then zoom in or out with a pinch or pull. Touch  (Back) to navigate to the *Post-acquisition screen*. Or, touch  (Trash) and then touch **Delete** at the dialogue screen: *Do you wish to discard this image and acquire again?* to confirm the discard and to return to the *Acquisition screen* to acquire again. Touch  (Add image) to navigate to the *Acquisition screen* and take an additional image.

Note The initial captured image is saved.


Touch  (Discard) then touch **Yes** at the dialogue: *Do you wish to discard this image?* to discard the image. Alternatively, touch **No** at the dialogue screen to return to the *Post-acquisition screen* to save the initial captured image.


On the *Post-acquisition screen*, touch  (Trash) and then touch **Delete** at the dialogue screen: *Do you wish to discard this image and acquire again?* to confirm the discard and to return to the *Acquisition screen* to acquire again.


Touch  (Accept) to accept the image and to navigate to the next image. Or touch  (Accept) to navigate to the *Exam summary screen* once the last image of the exam is acquired.


Exam summary controls (Post-acquisition and Review modes)


Exam summary from the *Post-acquisition* screen

From the *Post-acquisition* screen, touch  (Accept) to navigate to the *Exam summary* screen. Each image on the *Exam summary* screen contains a quality score and status. **(Wireless workflow images shown)**

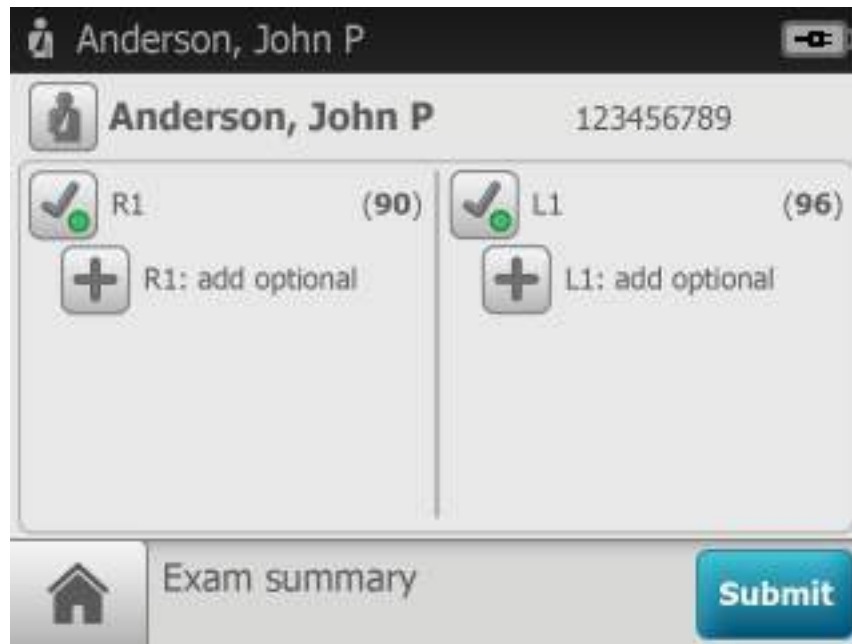

Touch  (Good quality green) to view that image. Image accepted, good quality. Quality scores of 40+.

Touch  (Acceptable quality yellow) to view that image. Image accepted, quality OK. Quality scores ranging from 20 – 39.


Touch  (Unacceptable quality red) to view that image. Image rejected. Quality scores ranging from 0 – 19.




CAUTION Captured images require patient context before the images can be saved.

Note

To add missing patient context information touch .

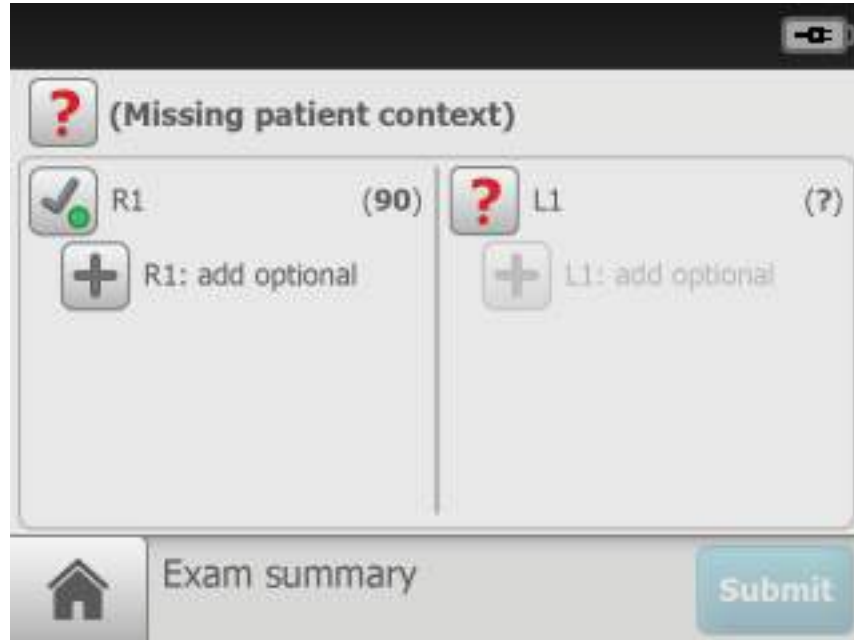

Touch  to create a new patient, or touch on an existing patient's name from within the *Patients* screen to select that patient.


(USB only) Touch **SAVE** to accept all images and to return to the *Home* screen.

(Wireless only) Touch **Submit** to accept and send all images and to return to the *Home* screen.

If an image is missing touch **Override** to access the Submit (or Save) button and to submit or save a single image in the 1-image, Macula-centered Exam setup,




or touch **Override** to submit (or save) multiple images in the 2-image, Macula + Optic Disc-centered Exam setup.

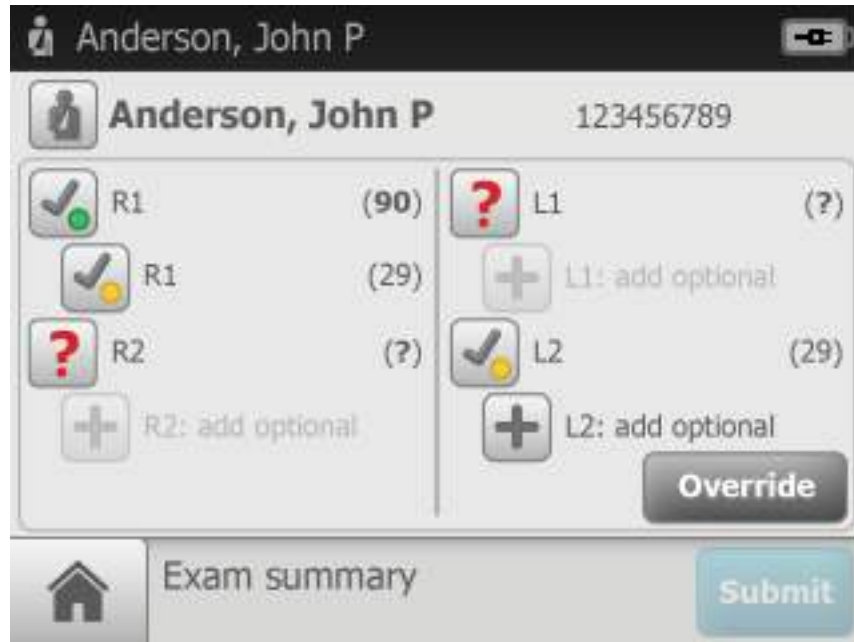





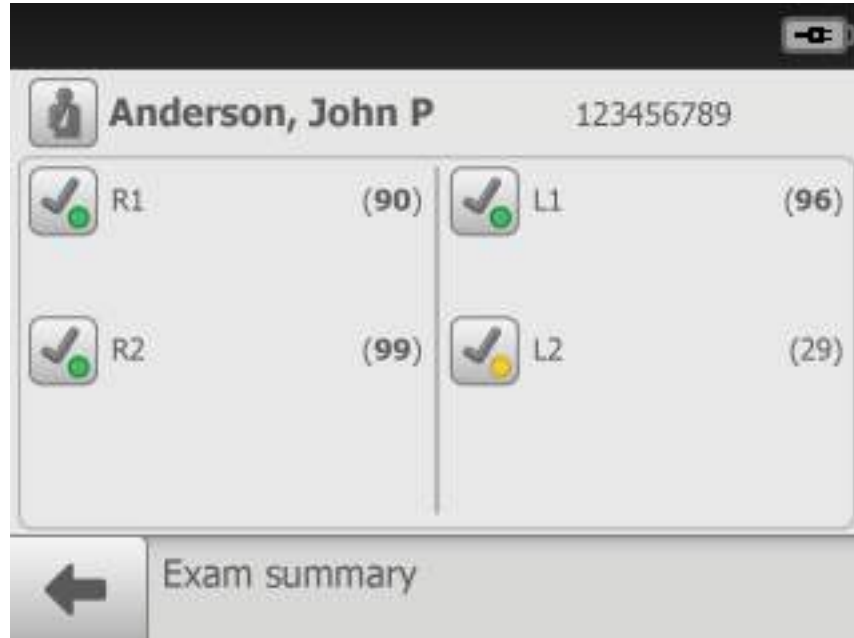

Exam summary from the *Review screen*

Touch **Review**.

Touch on a patient's name. The *Exam summary screen* appears for that patient.




Touch  (Good quality green),  (Acceptable quality yellow), or  (Unacceptable quality red) to view the image. Touch  (Patient) to review patient

demographics on the *Patient summary* screen. Touch  (Back) to return to the *Review screen*. Touch  (Home) to return to the *Home screen*.

Power down the RetinaVue 100 Imager

1. Press .

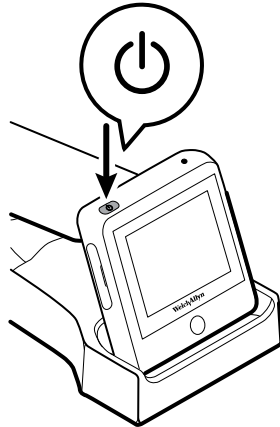

2. A dialogue box appears with three options:
 - Power Off
 - Sleep
 - Cancel




3. Touch one of the options.
The RetinaVue 100 Imager will power off, go into Sleep mode, or return to the prior screen, depending on the option you choose.

Sleep mode

To conserve power, the RetinaVue 100 Imager enters Sleep mode if approximately 2 minutes have passed since the last screen press.

To bring the RetinaVue 100 Imager out of Sleep mode, use any of these methods:


1. Touch the power button.
2. Touch the LCD touchscreen.
3. Press the Home button.

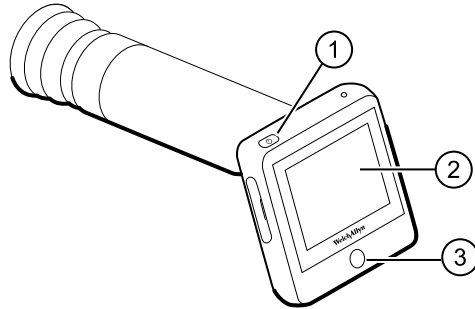

Note Lifting or moving the RetinaVue 100 Imager also brings the RetinaVue 100 Imager out of Sleep mode.

Enter Sleep mode

1. Press .

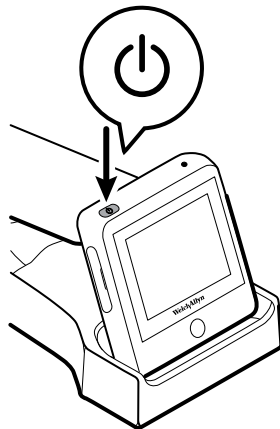

2. A dialogue box appears with three options:
 - Power Off
 - Sleep
 - Cancel




3. Touch **Sleep**.

The RetinaVue 100 Imager enters Sleep mode.

Exit sleep mode

1. Perform one of these options: (1) Press , (2) touch the LCD touchscreen, or (3) press the Home button.

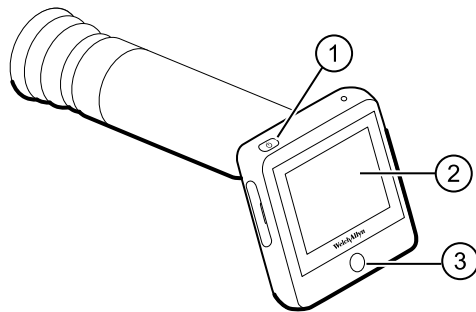

Note Lifting or moving the RetinaVue 100 Imager also brings the RetinaVue 100 Imager out of Sleep mode.

2. After performing one of the options, the RetinaVue 100 Imager will awaken from Sleep mode, or return to the prior screen.

Update the RetinaVue 100 Imager software to version 6.XX and above (USB and wireless)

Note (USB only) Patient list changes will only be made to the RetinaVue 100 Imager after it is undocked from the docking station charger. (Transferring patient list information requires additional RetinaVue Network software.)

Note (USB only) DICOM files will be available to the RetinaVue Network once the RetinaVue 100 Imager is docked.




CAUTION Do NOT remove the SD Card or the power source from the device during boot up.

Take note of the current software version:

1. Touch **Settings**.
2. Touch **About** to view the software version.

Before you begin, you will need the following :

- A computer with an Internet connection.
 - The docking station charger with a USB cable connected to the power source.
 - The RetinaVue 100 Imager containing the SD card that you plan to update. (Alternatively, you can use a blank, formatted, 8GB to 32GB SD Card to transfer the downloaded software update to the RetinaVue 100 Imager).
1. Ensure that any patient data that you want to save is transferred to the RetinaVue Network by docking the RetinaVue 100 Imager into the docking station charger with a USB cable connected to a laptop or PC. (If the RetinaVue Network software is not already open, double-click the desktop icon.)
 2. Remove the RetinaVue 100 Imager from the docking station charger, touch **Settings** then touch **Advanced settings**.
 3. Touch **Format SD**.
 4. Touch **OK**.
 5. Power down the RetinaVue 100 Imager and remove the SD Card.

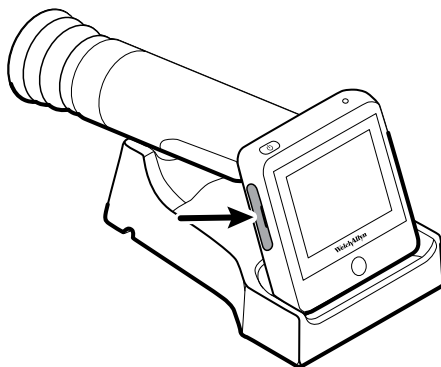

6. Insert the SD Card into an open SD Card reader port on the computer.
7. Download the software update from https://www.retinavue.net/RN_CustomerPortal/ and transfer the software update to a computer.
8. Transfer the software update from the computer to the SD Card.

Note Use a safe eject method when ejecting the SD card to ensure a complete software transfer to the SD card. (Options include right-clicking on the *Start Menu > Open Windows Explorer> right-click on SD card> click Eject*. Or, click on the *Taskbar> Show hidden icons > click Safely Remove Hardware and Eject Media> click Eject SDHC*). If an eject menu is not available with your system settings, wait until the files have completely copied before manually ejecting the SD Card.

9. Remove the SD Card from the computer.
10. Insert the SD Card into the RetinaVue 100 Imager.
11. Place the RetinaVue 100 Imager on the docking station charger, connect the docking station charger via a USB cable to the power source, and then power up the RetinaVue 100 Imager.

Note The update process typically takes a few minutes to cycle through the complete update. During the update process, it is normal to see both a blank screen and an update screen that appear several times before the RetinaVue 100 Imager automatically reboots.

12. After the RetinaVue 100 Imager boots up to the *Home* screen, confirm the software update.
 - a. Touch **Settings**.
 - b. Touch **About** to view the software version.
13. Touch **OK**.

Update the RetinaVue 100 Imager connection settings (USB and wireless)

Complete the following steps to update or change USB or wireless connection settings:

14. Touch **Advanced settings**.
15. Touch **Restore factory defaults**.
16. Touch **OK** at the prompt: *Reset all the settings to defaults*.
17. Touch **OK** at the prompt: *All settings were restored to factory defaults*.

You are ready to use the updated software. See either topic *First time selection of USB connectivity workflow* or *First time selection of Wi-Fi connectivity workflow* for further information about first time set up of USB or wireless connectivity.

Troubleshooting software updates

Note The RetinaVue 100 Imager boots up normally after a successful software update, however if a software update fails, the RetinaVue 100 Imager displays a failed update error message.

Power loss during software update:

- If a power loss occurs during the software update or if there are hardware issues, the RetinaVue 100 Imager posts a blank screen. Remove the battery, reinstall it, and restart RetinaVue 100 Imager. By performing this power cycle, the RetinaVue 100 Imager reverts to a previous version of software. Repeat the steps to perform the update again.

Corrupted software files:

- If the software update on the SD Card is corrupted, the RetinaVue 100 Imager displays a failed update error message indicating that the RetinaVue 100 Imager is unable to load the software from the SD Card. Remove the SD Card, reformat it on a PC. Repeat the steps to perform the update again.

Update extraction file error

If the RetinaVue 100 Imager displays the failed update error message, *Update file extraction failed. Power off device and restart*, the SD Card does not have sufficient space. Remove all files (except the update image) from the SD Card. Perform the update again.


Viewing or changing the RetinaVue 100 Imager settings


To view or change the settings




Touch **Settings**. Modify the settings as desired.

- Touch **WiFi On/Off**.
 - Touch (On) or (Off) to turn the wireless radio on or off.
- Touch **WiFi**. The list of available wireless networks in range appears.



Touch **>** to view or modify the wireless network settings from the list of available


wireless networks in range. Touch **OK** to save any changes. Touch  (Back) to return to the *Settings screen*.

- Touch **Clinic / States**. Up to 6 States may be selected for a clinic.
 - Touch **Multi-State**. The switch changes from Off to On.
 - Touch **States**. Touch the box next to the State abbreviation to select it.
 - Touch **More**. Touch **OK**.
 - Touch  (Back) to return to the *Settings screen*.
- Touch **Volume**.




Touch  (up) or  (down) to adjust the device sound volume. Touch  (Back) to return to the *Settings screen*.

- Touch **Time**.

Touch  (up) or  (down) to adjust the hour, minute, Continent, and Location.

Use the radio buttons to select the 24-hour or 12-hour format. Touch  (Back) to return to the *Settings screen*.


- Touch **Date**.




Touch  (up) or  (down) to adjust the day, month, and year in the DD/MM/YY format. Touch  (Back) to return to the *Settings screen*.


Touch **Advanced settings**.

Modify the settings as desired

- **Device access code**. Set the Device access numeric code, retype the code to confirm, and touch **OK** to return to the *Advanced settings screen*. (USB only)
- **Advanced settings code**. Set the Advanced settings numeric code, retype the code to confirm, and touch **OK** to return to the *Advanced settings screen*.
- **Format SD**. Touch **OK** to format the SD card to the default settings and touch **OK** again to return to the *Advanced settings screen*.
- **Restore factory defaults**. Touch **OK** to reset the device to the default settings and touch **OK** to return to the *Advanced settings screen*.
- **Exam setup**. Touch **Exam setup**. Touch 1-image (Macula-centered) to use the one image per eye exam workflow or touch 2-image (Macula and Optic Disc-centered) to

use the two images per eye exam workflow . Touch  (Back) to return to the *Settings* screen.


- **Export Lg.** The log files are exported to the SD card. Touch **OK** to return to the Advanced Settings screen. Touch  (Back) to return to the *Settings* screen.
- Touch **About** to view the software version. Touch **OK**.
- **Regulatory Information.** Touch **Regulatory Information** to view Regulatory Information. Touch  (Back) to return to the *Settings* screen.
- **Network Information.** Touch **Network Information** to view the wireless radio (Wi-Fi) MAC address, IP address, Netmask address, or Broadcast address. Touch  (Back) to return to the *Settings* screen.

Touch  (Home) to return to the *Home* screen.


Set the date and time

Before first use, set the date and time:

1. Follow the steps in the *Battery installation and replacement* section to install the battery and then press and hold the power button to turn the device on.

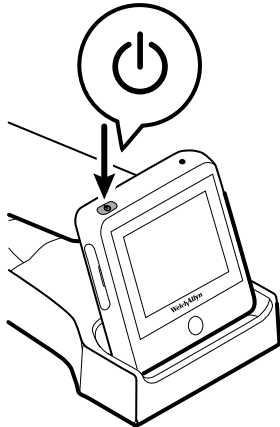




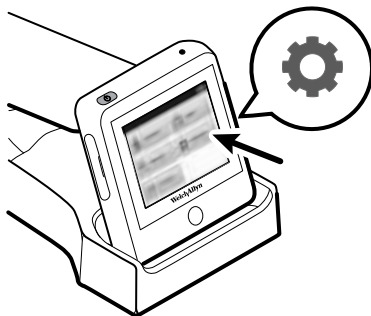


2. Touch **Settings** on the RetinaVue 100 Imager Home screen.

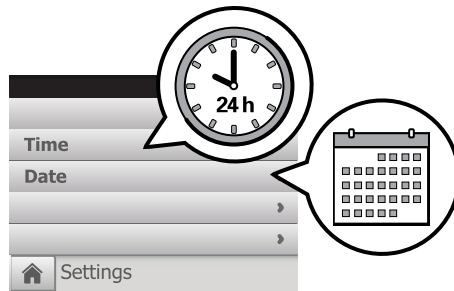

3. Touch **Time**. Touch  (up) or  (down) to adjust the hour, minute, Continent, and location.

Note The default time appears in the 24-hour format. Use the radio buttons to select the 24-hour or 12-hour format.




4. Touch  (Back) to return to the *Settings* screen.
5. Touch **Date**. Touch  (up) or  (down) to adjust the day, month, and year in the DD/MM/YY format.
6. Touch  (Back) to return to the *Settings* screen.
7. Touch  (Home) to return to the RetinaVue 100 Imager *Home* screen.

Set the Device access code at the RetinaVue 100 Imager camera (USB only)


Note Setting the Device access code before first use provides additional security features in the RetinaVue 100 Imager to prevent unintended access to patient information by unauthorized individuals.

Note See *First time RetinaVue 100 Imager camera registration at the RetinaVue Network Customer Portal* for instructions on setting the Device access code on the camera from the RetinaVue Network Customer Portal.

1. Follow the steps in the *Battery installation and replacement* section to install the battery and then press and hold the power button to turn the device on.

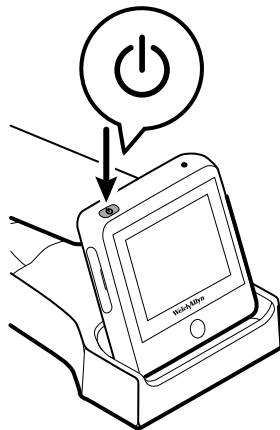

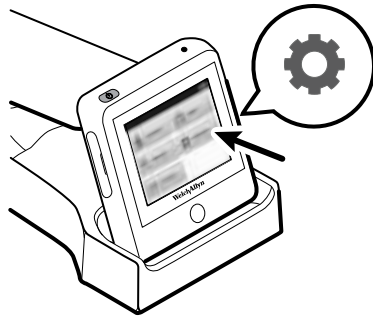

2. Remove the RetinaVue 100 Imager from the docking station charger.
3. Touch **Settings** on the RetinaVue 100 Imager Home screen.



4. Touch **Advanced settings**.
5. Touch **Device access code** to enter an access code for the RetinaVue 100 Imager.
6. Use the touchscreen numeric keypad to enter an access code into the *Access code* field and then retype the same code into *Access code (retype)* field.


Note The access code password can contain a minimum of 1-digit and a maximum of 18-digits. Characters are only numeric from 0 through 9.

7. Touch **OK** to accept the code and touch **OK** again at the confirmation screen to return to the *Advanced settings* screen.
8. Touch  (Back) to return to the *Settings screen*.
9. Touch  (Home) to return to the RetinaVue 100 Imager *Home screen*.

Reset the Device access code (USB only)




Note See the *Troubleshooting section* if you cannot recall your RetinaVue 100 Imager access code and need to erase the access code.

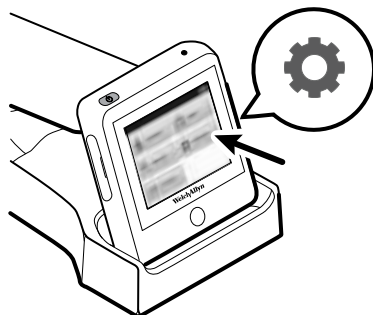1. Remove the RetinaVue 100 Imager from the docking station charger.
2. If you are not already logged in, enter the existing access code for the RetinaVue 100 Imager.
3. Touch **Settings** on the RetinaVue 100 Imager Home screen.





4. Touch **Advanced settings**.
5. Touch **Device access code**.
6. Touch **Reset** to clear the existing access codes.
7. Use the touchscreen numeric keypad to enter an access code into the *Access code* field and then retype the same code into the *Access code (retype)* field.

Note The access code password can contain a minimum of 1-digit and a maximum of 18-digits. Characters are only numeric from 0 through 9.

Note Alternatively, after the code reset, touch  (Back),  (Home), and  (Back) to return to the RetinaVue 100 Imager *Home screen* without using a Device access code.

8. Touch **OK** to accept the code and touch **OK** again at the confirmation screen to return to the *Advanced settings* screen.


9. Touch  (Back) to return to the *Settings screen*.

10. Touch  (Home) to return to the RetinaVue 100 Imager *Home screen*.

Reset the Advanced settings code

Note See the *Troubleshooting section* if you cannot recall your RetinaVue 100 Imager Advanced settings code or the Device access code and need to erase the code.

1. Remove the RetinaVue 100 Imager from the docking station charger.
2. If the RetinaVue 100 Imager is set with a Device access code (in addition to the Advanced settings code) and you are not already logged in, enter the existing Device access code.
3. Touch **Settings** on the RetinaVue 100 Imager *Home screen*.

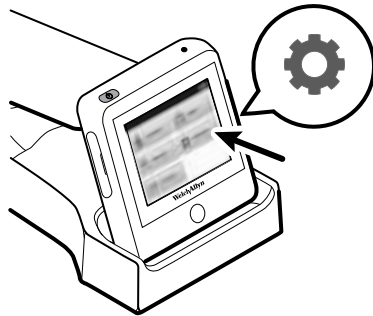


4. Touch **Advanced settings**.
5. Enter the existing Advanced settings code for the RetinaVue 100 Imager.
6. Touch **Advanced settings code**.
7. Touch **Reset** to clear the existing Advanced settings codes.
8. Use the touchscreen numeric keypad to enter a code into the *Advanced settings code* field and then retype the same code into the *Advanced settings code (retype)* field.


Note The Advanced settings code password can contain a minimum of 1-digit and a maximum of 18-digits. Characters are only numeric from 0 through 9.

Note


Alternatively, after the code reset, touch

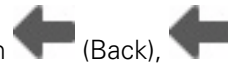


(Back),

 (Home) to return to the RetinaVue 100 Imager *Home screen* without using an Advanced setting code.

9. Touch **OK** to accept the code and touch **OK** again at the confirmation screen.

10. Touch  (Back) to return to the *Settings screen*.

11. Touch  (Home) to return to the RetinaVue 100 Imager *Home screen*.

RetinaVue 100 Imager Accessories

Miscellaneous items

Part number	Description
RV100-B	RetinaVue 100 Imager with lithium-ion battery
RV100-6	RetinaVue 100 Imager with lithium-ion battery and Australian plug
106401	Charging Dock for Welch Allyn RetinaVue 100 Imager; 100-240 V, 50-60 Hz
106405	Battery for Welch Allyn RetinaVue 100 Imager; Lithium-Ion 3.6V / 3350mAh
106406	Replacement USB Cable for Welch Allyn RetinaVue 100 Imager
106507	Power supply wall adapter for Welch Allyn RetinaVue 100 Imager. 100-240 V, 50-60 Hz; IEC Plug Type B (US)
106508	Power supply wall adapter for Welch Allyn RetinaVue 100 Imager. 100-240 V, 50-60 Hz; IEC Plug Type C (EU)
106509	Power supply wall adapter for Welch Allyn RetinaVue 100 Imager. 100-240 V, 50-60 Hz; IEC Plug Type G (UK)
106530	Power supply wall adapter for Welch Allyn RetinaVue 100 Imager. 100-240 V, 50-60 Hz; IEC Plug Type I (AU)
106710	RetinaVue 100 Imager Lens cap
106711	RetinaVue 100 Imager Battery cover
236200-2	Welch Allyn RetinaVue 100 Imager replacement screwdriver for battery cover. #1 Phillips Screwdriver w/ Logo
723481	Laptop for RetinaVue Network
724224	RetinaVue 100 Imager camera case
723484	Lens spray
723671	Lens cloth

Applied parts

Part number	Description
106404	Reusable Eye Cups for Welch Allyn RetinaVue 100 Imager: Qty. 5

Service

Part number	Description
S1-RV100-1	Welch Allyn Partners in Care Comprehensive Partnership Program Service for RetinaVue 100 Imager; 1 year
SV-RV100-2	Welch Allyn Partners in Care Comprehensive Partnership Program Service for RetinaVue 100 Imager; 2 years
SV-RV100-3	Welch Allyn Partners in Care Comprehensive Partnership Program Service for RetinaVue 100 Imager; 3 years
S1-RV100-1P	Welch Allyn Partners in Care Comprehensive Partnership Program Premium Service for RetinaVue 100 Imager; 1 year
S1-RV100-2P	Welch Allyn Partners in Care Comprehensive Partnership Program Premium Service for RetinaVue 100 Imager; 2 years
S1-RV100-3P	Welch Allyn Partners in Care Comprehensive Partnership Program Premium Service for RetinaVue 100 Imager; 3 years

Literature/Documentation

Part number	Description
724146	RetinaVue 100 Imager Technique guide
726974	RetinaVue Network Exam Wall Poster

Warranty

Welch Allyn warrants the product to be free of defects in material and workmanship and to perform in accordance with manufacturer's specifications for the period of one year from the date of purchase from Welch Allyn or its authorized distributors or agents.

The warranty period shall start on the date of purchase. The date of purchase is: 1) the invoiced ship date if the device was purchased directly from Welch Allyn, 2) the date specified during product registration, 3) the date of purchase of the product from a Welch Allyn authorized distributor as documented from a receipt from said distributor.

This warranty does not cover damage caused by: 1) handling during shipping, 2) use or maintenance contrary to labeled instructions, 3) alteration or repair by anyone not authorized by Welch Allyn, and 4) accidents.

Shipping cost to return a device to a Welch Allyn Service center is not included.

A service notification number must be obtained from Welch Allyn prior to returning any products or accessories to Welch Allyn's designated service centers for repair. To obtain a service notification number, contact Welch Allyn Technical Support.

THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. WELCH ALLYN'S OBLIGATION UNDER THIS WARRANTY IS LIMITED TO REPAIR OR REPLACEMENT OF PRODUCTS CONTAINING A DEFECT. WELCH ALLYN IS NOT RESPONSIBLE FOR ANY INDIRECT OR CONSEQUENTIAL DAMAGES RESULTING FROM A PRODUCT DEFECT COVERED BY THE WARRANTY.

WelchAllyn[®]

Advancing Frontline Care[™]

Material No.

411492